


# Syringe Guide



The Measure of Excellence®  
for your application

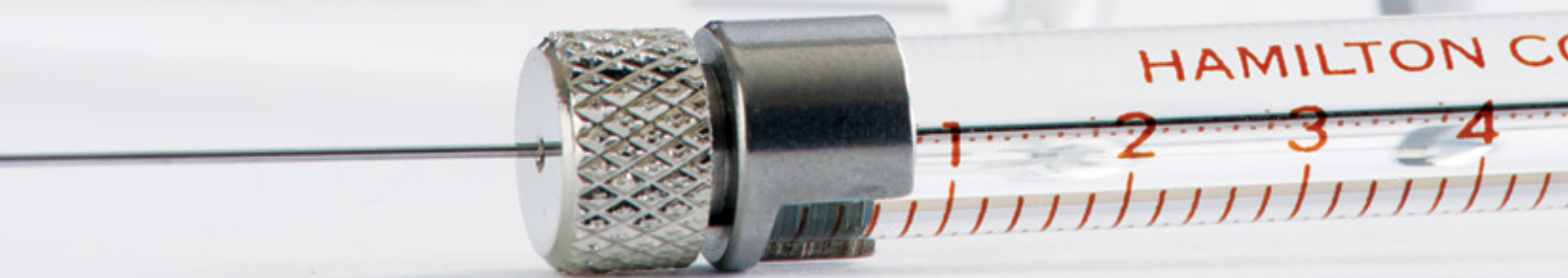

**HAMILTON** 

# Unsurpassed Accuracy and Precision

Hamilton syringes are the finest quality precision fluid measuring devices available. Top quality materials and skilled workmanship ensure that Hamilton syringes consistently deliver the highest possible performance for reliable analyses. For manual uses, our syringes are accurate to within  $\pm 1\%$  of the nominal volume with a precision of 1% at 80% of the total volume.

## A New Application-Based Guide

Hamilton offers the most complete selection of syringes on the market for use in various applications including gas chromatography (GC) and high performance liquid chromatography (HPLC), liquid handling and life sciences. The present guide is a tool to help you select the appropriate syringe for your application. It also provides detailed background and technical information to help you understand the full functionalities and possibilities for your syringes.


**Products manufactured by Hamilton Company are intended for scientific research and laboratory use only and are not intended for human *in vivo* use.**



## Find the Perfect Syringe Online

If you do not find the correct syringe to meet your needs in this guide, please consult Hamilton's online syringe selection tool at [www.hamiltoncompany.com/syringes](http://www.hamiltoncompany.com/syringes) or contact customer support using the information on the back of the guide.



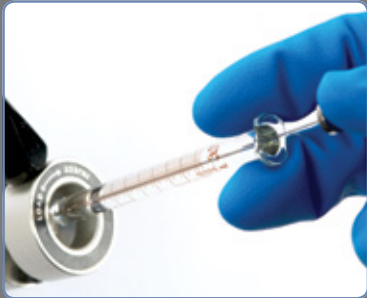

## HPLC autosampler syringes

See more on [page 20](#)



## Neuroscience animal injections

See more on [page 30](#)



## Syringe series and family definitions

See more on [page 4](#)

# Table of Contents

## Technical Syringe Information.....2

- Syringe Technical Illustrations .....2
- Syringe Series .....4
- Syringe Terminations .....4
- Needle Point Styles .....4

## Gas Chromatography Syringes .....5

- Manual Split/Splitless Injection .....6
- Manual On-Column Injection ..... 10
- Gas Chromatography Autosampler Syringes..... 11

## Liquid Chromatography Syringes ..... 16

- HPLC Syringes ..... 16
- Thin Layer Chromatography ..... 19
- HPLC Autosampler Syringes .....20

## Sample Preparation and Liquid Handling .....26

- Syringes For the Hamilton Microlab 600.....27
- Syringes For Laboratory Syringe Pumps .....27
- Tubing Assemblies .....27

## Life Science Syringes .....29

- General Animal Injections .....29
- Neuroscience Injections .....30
- Genomics and Biochemistry .....34

## Product Listing by Features.....36

- Microliter Syringes .....36
- Gastight Syringes .....38
- Needles.....40

## Syringe Care and Use.....44

## Part Number Index.....47

## About Hamilton Company .....49


For more information on the full portfolio of Hamilton syringes, needles and accessories or to order a product, please visit [www.hamiltoncompany.com](http://www.hamiltoncompany.com) or refer to the back of this guide for additional contact details.



# Technical Syringe Information


## Microliter™ Syringes (700 Series)

Microliter syringes are used with liquid only. They feature a stainless steel plunger which is individually fitted to the glass barrel. The result is a liquid tight seal with nearly frictionless movement. These syringes are ideal for organic samples that are not prone to precipitation, crystallization or bonding with glass.

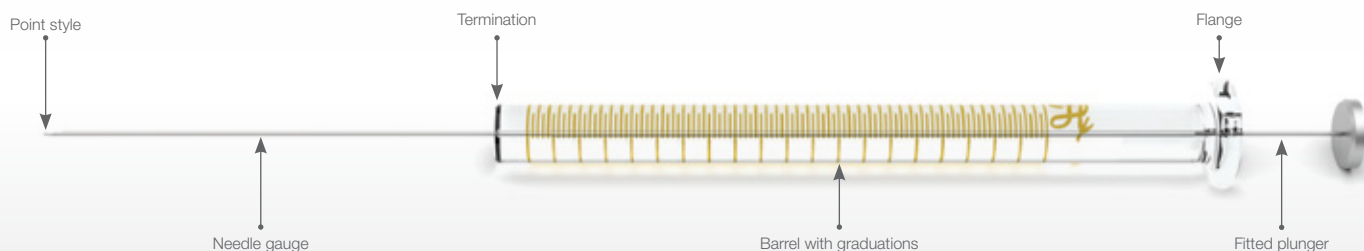


## Modified Microliter Syringes (7000 Series)

These syringes are designed to dispense volumes below 5  $\mu\text{L}$ . This requires an extremely small tungsten plunger wire that travels to the tip of the needle. The resulting design has zero dead volume and is accurate below 50 nL.

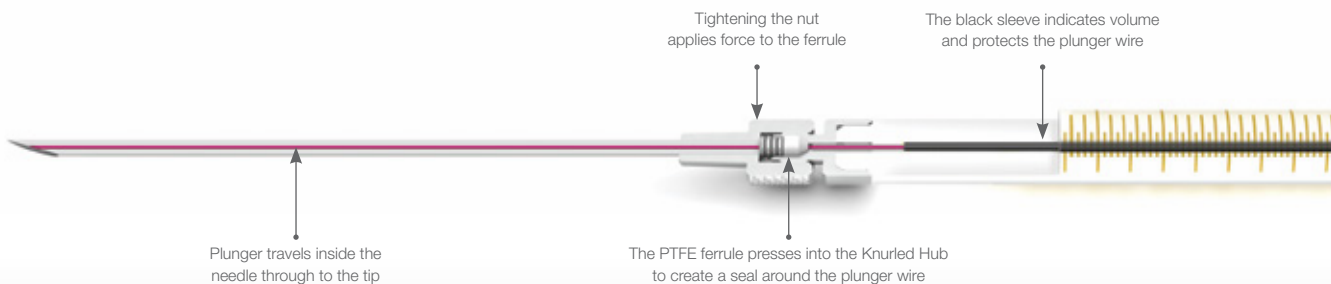


### Gastight® Syringes (1700 and 1000 Series)

Gastight syringes are ideal for dispensing both liquids and gases. They have a polymer plunger tip which creates a leak-free seal. Traditionally the tip is made from PTFE but other materials are used for selected applications. The polymer tip essentially wipes the interior of the syringe barrel free of sample. This feature is useful for aqueous and low volatile organic samples because it reduces the chance a deposit will occur and result in cross-contamination or even a damaged plunger.

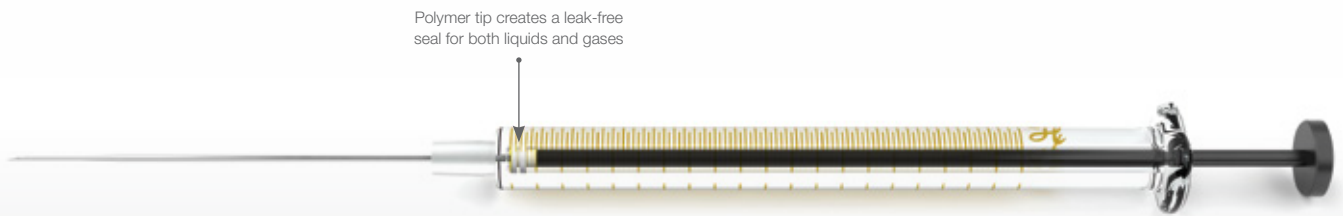

### Reinforced Plunger Syringes (800 and 1800 Series)




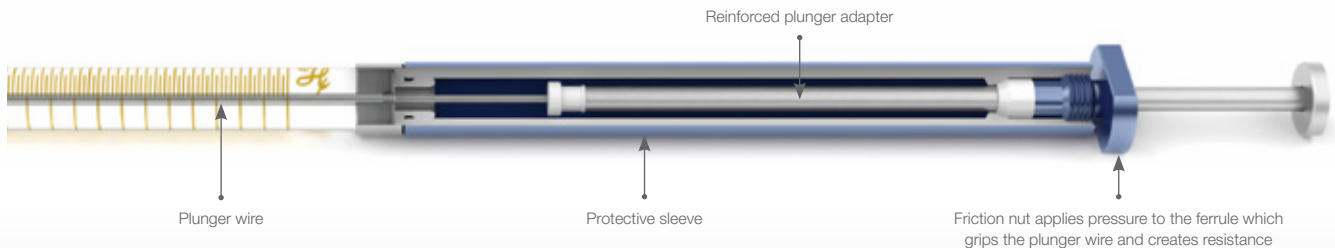





Small volume syringes require extremely small plunger wires. For uses that are prone to bending and breaking these small plungers, it is recommended to use a reinforced plunger syringe. These plungers are available in both Microliter and Gastight versions.

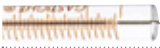







## Syringe Series

Syringe Series	Volume	Syringe Type	Plunger Type
7000	0.5 µL – 5 µL	Modified Microliter	Tungsten plunger
700	5 µL – 500 µL	Microliter	Stainless steel fitted plunger
800	5 µL – 250 µL	Microliter	Reinforced stainless steel fitted plunger
1700	10 µL – 500 µL	Gastight	Polymer tipped plunger
1800	10 µL – 250 µL	Gastight	Reinforced polymer tipped plunger
1000	1 mL – 100 mL	Gastight	Polymer tipped plunger
Super Syringe	0.5 L – 2 L	Gastight	Plunger with polymer tip with Buna-N seal

## Syringe Terminations

Termination	Image	Abbreviation	Remark	Autoclavable
Cemented Needle		N	For low-volume syringes	No
Luer Tip Cemented Needle		LTN	For mid-volume syringes	No
Removable Needle		RN		Yes if disassembled
Knurled Hub		KH	For Modified Microliter syringes only	Yes if disassembled
Luer Tip		LT		Yes if disassembled
PTFE Luer Lock		TLL		Yes if disassembled
SampleLock™		SL	Integrated on/off valve	No
Fixed Needle		FN	Found on CTC C-Line and X-Type syringes	No
ChemSeal™		C	1/4"-28 UNF/threaded connection	No


## Needle Point Styles

Point Style	Image	Description	Application
2		10–12° sharp, beveled, curved non-coring	Gas chromatography, septum piercing
3		Blunt, electro-polished	High performance liquid chromatography (HPLC) injection, thin layer chromatography (TLC), general liquid handling, controlled animal injections
3T		Blunt, electro-polished, coated with PTFE 19 mm from the tip	Thin layer chromatography (TLC) applications
4		Sharp 10–12° beveled needle	Life science/animal injections
5		Conical with side port for penetration without coring	Headspace, applications prone to needle clogging, causes minimal septum damage
AS		Conical, non-coring designed to withstand multiple injections	Autosampler injection, pre-pierced septa



# Gas Chromatography Syringes

The invention of the Microliter syringe by Clark Hamilton in 1947 enabled the development of gas chromatography. “At that time there was no convenient way to introduce microliter (or even smaller) volumes of liquids into a gas chromatograph,” wrote Professor Georges Guiochon, Professor and Distinguished Scientist at University of Tennessee. “In fact, it is not an exaggeration to say that without the Hamilton syringes, gas chromatography could not have become everybody’s technique” [Chromatographia 15(6), 333 (1982)]. Injection into a gas chromatograph may be manual or automatic. In both cases, a microsyringe is required to introduce the sample. The choice of the correct syringe depends on the type of injector used.

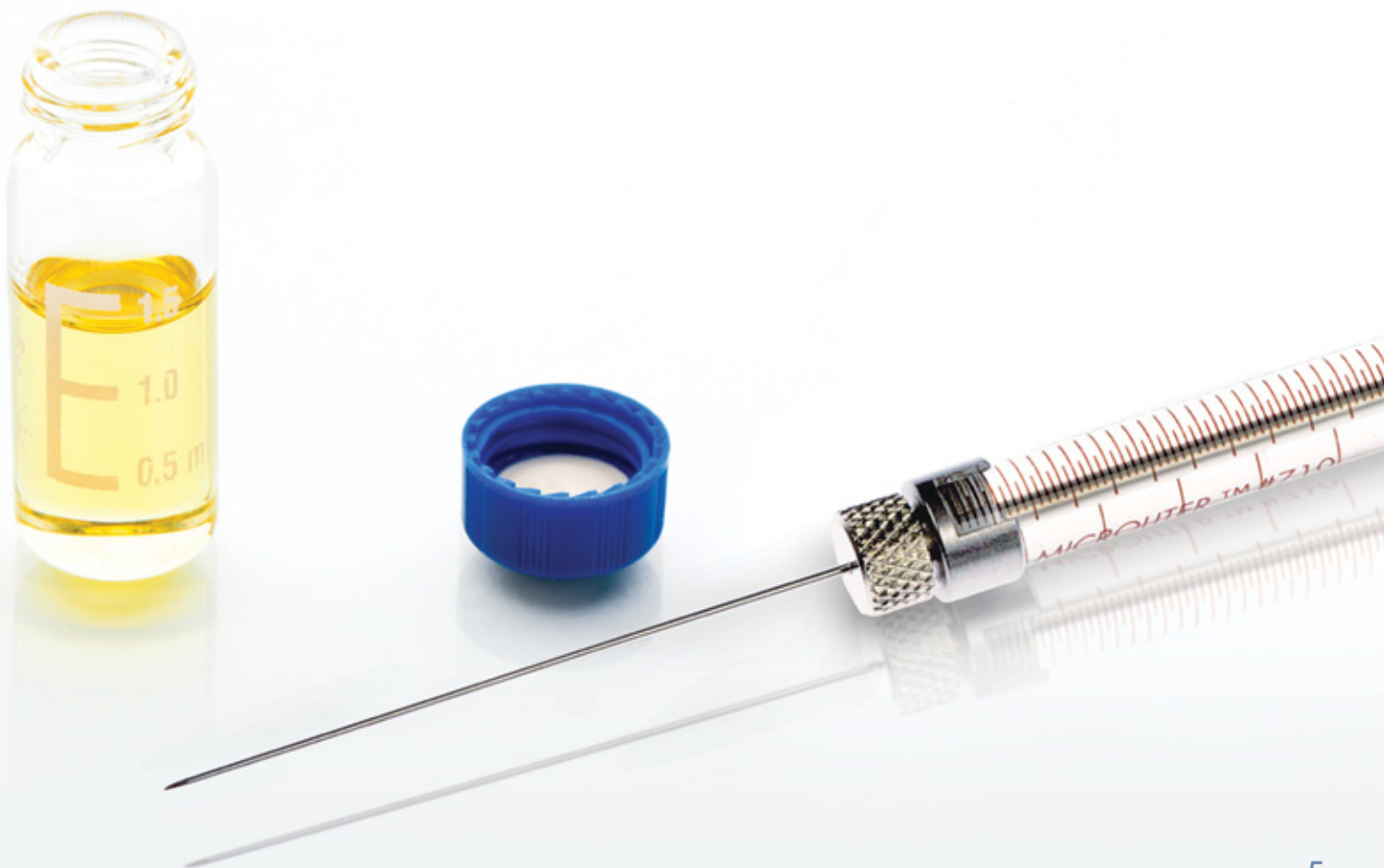

# Manual Split/Splitless Injection

A split/splitless injector consists of a heated chamber with a glass liner into which the sample is injected through the septum. A microsyringe is used to inject the sample through a rubber septum into a flash vaporizer chamber at the head of the column. The sample vaporizes to form a mixture of carrier gas, vaporized solvent and vaporized analytes. In the split mode, only a proportion of this mixture reaches the column but most exits through the split outlet. This avoids overloading the column. In the splitless mode, the split vent is closed so that the vaporized analyte passes onto the column. This mode is more sensitive and adequate for trace analysis.

## Common syringe features for split/splitless injection include:

- ▶ Microliter syringes for liquid samples and Gastight syringes for gas and liquid samples
- ▶ Common volume range between 0.5 to 50  $\mu\text{L}$
- ▶ Needle with point style 2 to pierce the rubber septum
- ▶ Needle length of 51 mm to reach the middle of the glass liner. Splitless injection of small volumes may be performed with a 70 mm needle to deposit the sample close to the column entrance.


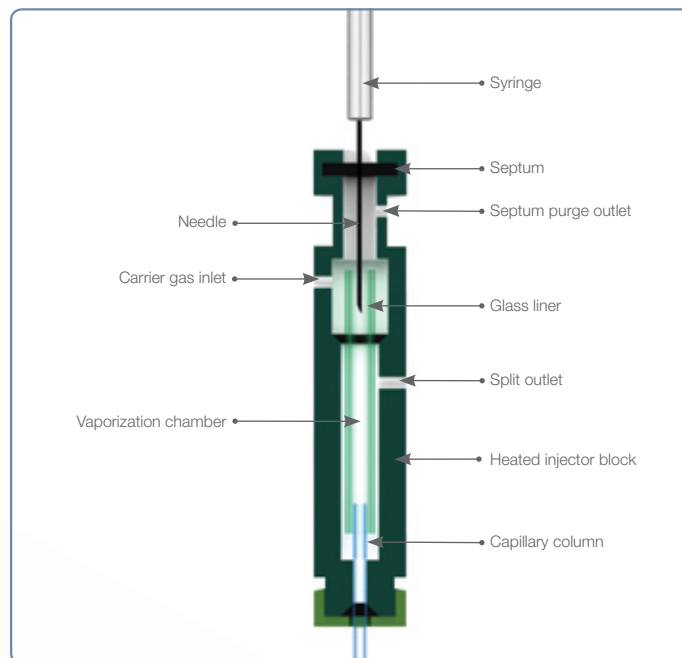


Illustration of a split/splitless injector



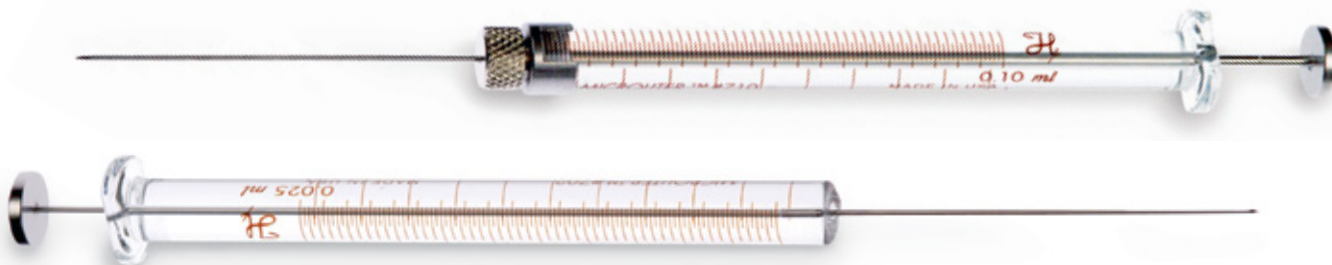

## Injection of Liquid Samples

### Syringes For Manual Injection of Liquid Samples


Volume	5 $\mu\text{L}$	10 $\mu\text{L}$	25 $\mu\text{L}$	50 $\mu\text{L}$	100 $\mu\text{L}$	250 $\mu\text{L}$	500 $\mu\text{L}$
Model	75	701	702	705	710	725	750
Gauge	26s	26s	22s	22s	22s	22s	22
Cemented Needle (N) Syringe (point style 2)	87900	80300 80366 (6 pk)	80400	80500	80600	80700	80800
Removable Needle (RN) Syringe (point style 2)	87930	80330 80336 (6 pk)	80430	80530	80630	80730	80830
Replacement Needle 6 pk (point style 2)	7758-02	7758-02	7758-03	7758-03	7758-03	7779-03	7779-01

Note: Needles are 51 mm

### Zero Dead Volume Syringes

Volume	0.5 $\mu\text{L}$	1.0 $\mu\text{L}$	2.0 $\mu\text{L}$	5.0 $\mu\text{L}$
Model	7000.5	7001	7002	7105
Gauge	25	25	25	24
Knurled Hub (KH) Syringe (point style 2)	86259	80135	88411	88011
Standard Injection Spacer (25 mm)	86201	86201	86201	86201

Note: Needles are 70 mm

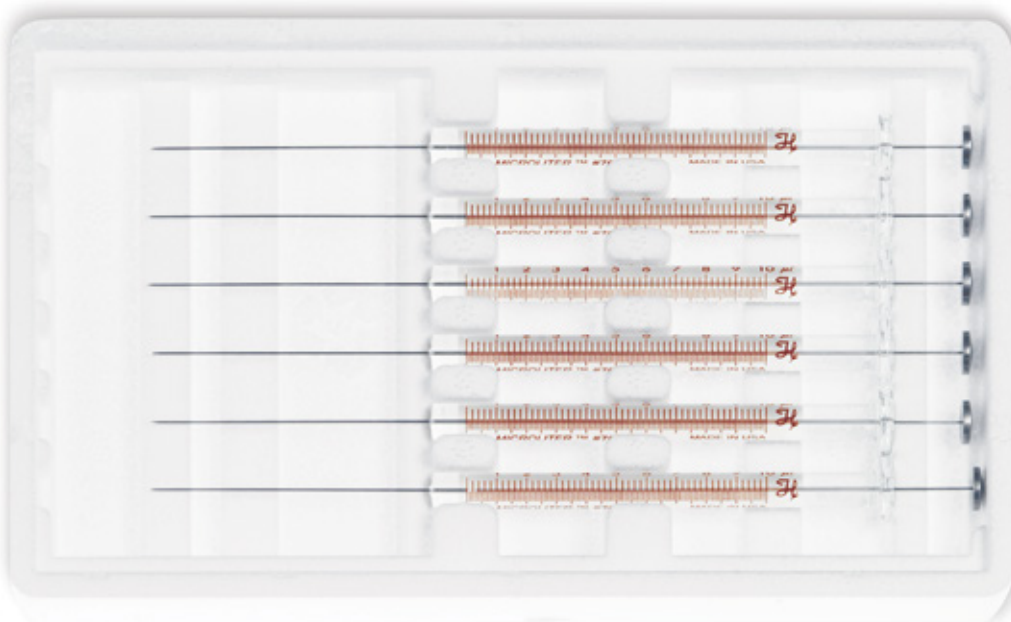

#### Tech Tip

Choose a syringe with a removable needle if working with halogenated solvents.

**Injection of Gas and Liquid Samples**

**Syringes For Manual Injection of Gas and Liquid Samples**

<b>Volume</b>	<b>10 µL</b>	<b>25 µL</b>	<b>50 µL</b>	<b>100 µL</b>	<b>250 µL</b>
<b>Model</b>	<b>1701</b>	<b>1702</b>	<b>1705</b>	<b>1710</b>	<b>1725</b>
<b>Gauge</b>	<b>26s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>
Cemented Needle (N) Syringe (point style 2)	80000	80200	80900	81000	
Cemented Needle (LTN) Syringe (point style 2)					81100
Removable Needle (RN) Syringe (point style 2)	80030	80230	80930	81030	81130
Replacement Needle 6 pk (point style 2)	7758-02	7758-03	7758-03	7758-03	7779-03
<b>Volume</b>	<b>500 µL</b>	<b>1000 µL</b>	<b>2500 µL</b>	<b>5000 µL</b>	<b>10000 µL</b>
<b>Model</b>	<b>1750</b>	<b>1001</b>	<b>1002</b>	<b>1005</b>	<b>1010</b>
<b>Gauge</b>	<b>22</b>	<b>22</b>	<b>22</b>	<b>22</b>	<b>22</b>
Cemented Needle (LTN) Syringe (point style 2)	81217	81317	81417	81517	81617
Removable Needle (RN) Syringe (point style 2)	81230	81330	81430	81530	81630
Replacement Needle 6 pk (point style 2)	7779-01	7779-01	7779-01	7779-01	7779-01

Note: Needles are 51 mm

**Reinforced Plunger Syringes**

<b>Volume</b>	<b>10 µL</b>	<b>25 µL</b>	<b>50 µL</b>	<b>100 µL</b>	<b>250 µL</b>
<b>Model</b>	<b>1801</b>	<b>1802</b>	<b>1805</b>	<b>1810</b>	<b>1825</b>
<b>Gauge</b>	<b>26s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>
Removable Needle (RN) Syringe (point style 2)	84877	84880	84883	84886	84889
Replacement Needle, 6 pk (point style 2)	7758-02	7758-03	7758-03	7758-03	7779-03

Note: Needles are 51 mm



Reinforced Plunger Syringe



### Adapter For Repeated Manual Injection

Volume	10 $\mu$ L	25 $\mu$ L	50 $\mu$ L	100 $\mu$ L	250 $\mu$ L	500 $\mu$ L
Model	701/1701	702/1702	705/1705	710/1710	725/1725	750/1750
Reproducibility (Chaney) Adapter	14700	14725	14725	14725	14725	14725



*Chaney Adapter*



### Tech Tip

Improve the reproducibility of multiple manual injections with the Chaney Adapter (p/n 14700 and 14725).



# Manual On-Column Injection

On-column injection employs the direct introduction of the liquid sample into the column. This is very useful for the analysis of thermally unstable samples and avoids boiling point discrimination. On-column injection requires special syringes with the following features:

- ▶ Small outer diameter needles adjusted to the inner diameter of the column (0.17 mm needle for columns with 0.25 mm inner diameter, 32 gauge for columns with 0.32 mm inner diameter, and 26 gauge for columns with a 0.53 mm inner diameter)
- ▶ Longer needles tailored to pass through the injector and reach within the capillary column
- ▶ Point style 3 adapted to the typical septum of the on-column injector
- ▶ Small volume of 5  $\mu\text{L}$  or 10  $\mu\text{L}$


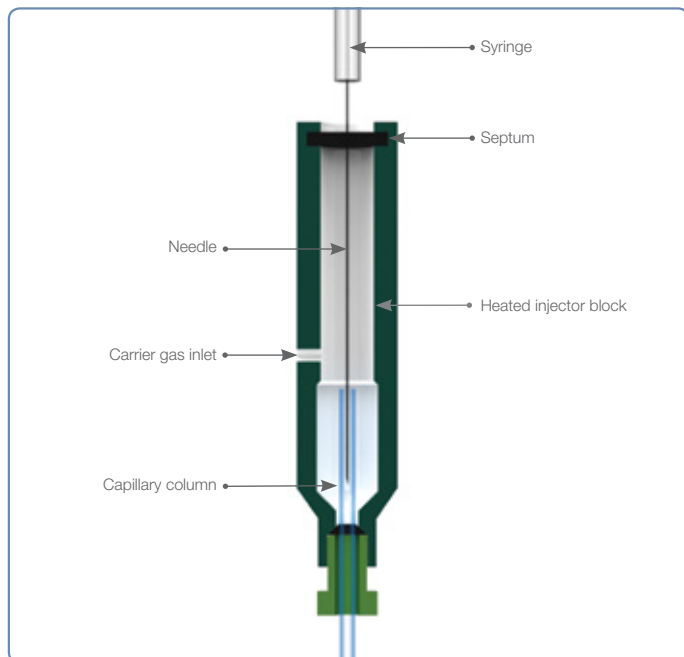


Illustration of an on-column injector

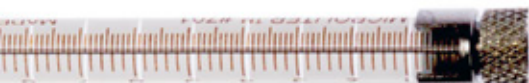

## Syringes for On-Column Injection For Capillary Columns

Volume	10 $\mu\text{L}$	10 $\mu\text{L}$	10 $\mu\text{L}$
Model	701	701	701
Gauge	0.17 mm	32	26s
Column I.D.	0.25 mm	0.32 mm	0.53 mm

Fused Silica Needle Syringe (point style 3)	87402 <sup>1</sup>	80308 <sup>2</sup>	80308 <sup>2</sup>
---	--------------------	--------------------	--------------------

Note: Needles are 100 mm

1) This is an RN needle

2) These are stainless steel needles. Please specify gauge, length and point style.




# Gas Chromatography Autosampler Syringes

Autosamplers enable the automatic introduction of samples into the injector of the gas chromatograph (GC). Automatic injection has become very common as it improves reproducibility and speed.

Hamilton offers a wide range of autosampler syringes. The features of the autosampler syringes are adapted for an exact fit to a specific autosampler. The point style AS (for autosampler) has been specifically designed to withstand repeated penetration through the GC autosampler septa. Syringes for headspace injection usually have a point style 5, which is ideal for large gas volume injection. Autosampler syringes are available in two gauge types depending on the type of injector used.

Gauge	Application
23s	Durable needle for Merlin MicroSeal™ septa and split/splitless injections
26s	Versatile needle for all types of injections including split/splitless and on-column injections

For a full list of all syringes available, visit [www.hamiltoncompany.com/syringes](http://www.hamiltoncompany.com/syringes)

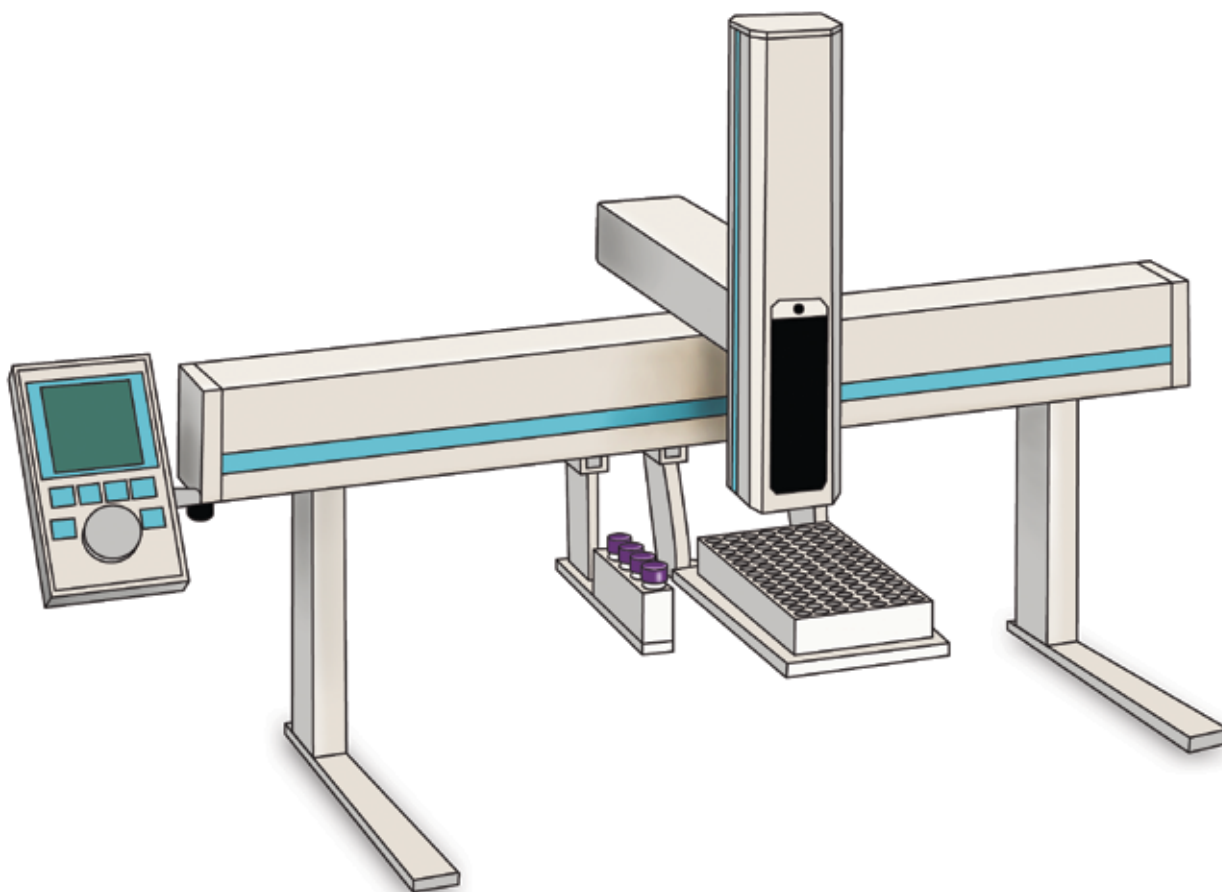

CTC Autosampler

## CTC PAL® COMBI-xt and PAL® GC-xt Autosampler Syringes

CTC PAL® autosamplers are one of the most popular sampling devices on the market today for GC and headspace analyses. Customers rely on its ability to sample from many different vial and well types and on its ultimate platform flexibility. The syringe is at the heart of every injection that the GC PAL® system performs and CTC chose Hamilton to provide this critical component.

### C-Line Syringes

C-Line syringes incorporate several unique design features that ensure superior performance. The innovative direct attachment design of the needle to the barrel eliminates contact between the sample solvent and adhesive, reducing carryover and ghost peaks. The inert plunger tip polymer provides enhanced solvent compatibility and longer lifetime.

### C-Line PAL® COMBI-xt and PAL® GC-xt Liquid Injection Syringes



Volume	1.2 µL	5 µL	10 µL	25 µL	100 µL	250 µL	500 µL
Model	7701.2	75	701	1702	1710	1725	1750
Gauge	26	26s	26s	26s	26s	26	26
Fixed Needle (FN) Syringe (point style AS)	203185 <sup>1</sup>	203189	203205 203361 <sup>2</sup>	203043	203076 <sup>3</sup>	203078	203080

Note: Needles are 51 mm

1) KH Termination

2) 23s gauge

3) Slimline syringe

### Headspace Syringes for CTC PAL® Combi-xt

Modern GC headspace analysis requires injecting over large temperature ranges. Conventional headspace syringes on the market use a rubber O-ring sealed plunger which has a limited sealing performance at high temperatures due to varying thermal expansion between the different materials. The high dynamic HD-type syringe employs a unique spring in the plunger tip which compensates for the materials' different expansion coefficients, creating a better seal over a larger temperature range, and it also improves lifetime.



### PAL® COMBI-xt Headspace Syringes

Volume	1000 µL	1000 µL	2500 µL	2500 µL	5000 µL	5000 µL
Model	1001	1001	1002	1002	1005	1005
Gauge	26	23	26	23	26	23
Cemented Needle (LTN) Syringe (point style 5)	203141	203082	203181	203084	203182	203086

Note: Needles are 51 mm



## Agilent 7673, 7683, 7693 and 6850 Autosampler Syringes

### Durable Syringes For Merlin Microseal™ Septa and Split/Splitless Injection


Volume	0.5 µL	1 µL	5 µL	10 µL	10 µL
Model	7000.5	7001	75	701	1701
Gauge	23	23	23s	23s	23s
Cemented Needle (N) Syringe (point style AS)			87987 87990 (6 pk)	80387 80390 (6 pk)	80080 80094 (6 pk)
Knurled Hub (KH) Syringe (point style AS)	86276	80176			

Note: Needles are 43 mm

### Universal Syringes For Split/Splitless and On-Column Injection

Volume	5 µL	10 µL	10 µL
Model	75	701	1701
Gauge	23s-26s	23s-26s	23s-26s
Cemented Needle (N) Syringe (point style AS)	87993 87994 (6 pk)	80393 80391 (6 pk)	80079 80096 (6 pk)


Note: Needles are 43 mm

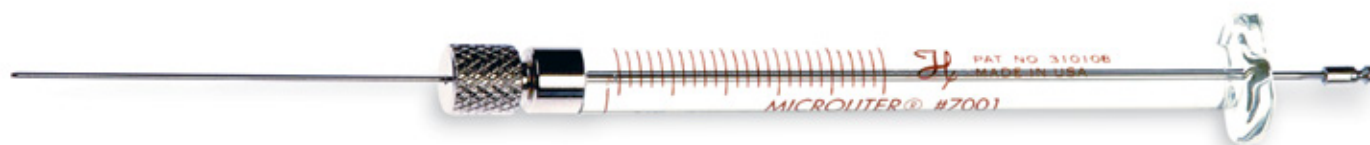

### Versatile Syringes For All Types of Injections, Including On-Column Injection

Volume	0.5 µL	5 µL	10 µL
Model	7000.5	75	701
Gauge	26	26s	26s
Cemented Needle (N) Syringe (point style AS)	86274 <sup>1</sup>	87988 87989 (6 pk)	80388 80389 (6 pk)

Note: Needles are 43 mm  
1) KH termination

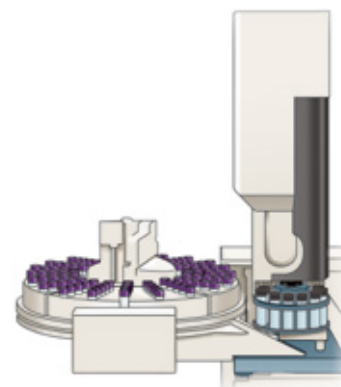

Agilent 7693 Autosampler

## Thermo Scientific® Autosampler Syringes

### AS 800, AS 2000 Syringes

Volume	10 µL
Model	701
Gauge	26s

Cemented Needle (N) Syringe (point style AS)	80318
--	-------


Note: Needles are 80 mm

### Headspace HS 250/500/850 Syringes

Volume	2500 µL
Model	1002
Gauge	22

Cemented Needle (LTN) Syringe (point style 5)	202660
---	--------

Note: Needles are 56 mm

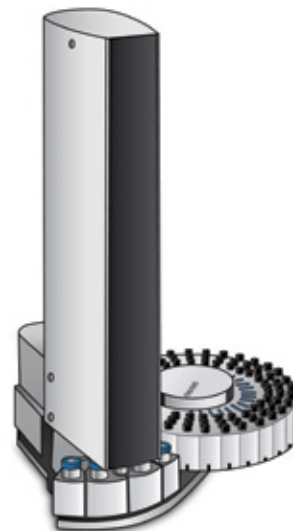

Thermo Scientific® Finnigan AS 2000 Autosampler

### AI/AS 3000 Tri Plus Syringes

Volume	5 µL	10 µL
Model	75	701
Gauge	26s	26s


Cemented Needle (N) Syringe (point style AS)	204051	204052
--	--------	--------

Cemented Needle (N) Syringe (point style 2)	204000	204001
---	--------	--------

Note: Needles are 50 mm



Varian 8100/8200 Autosampler

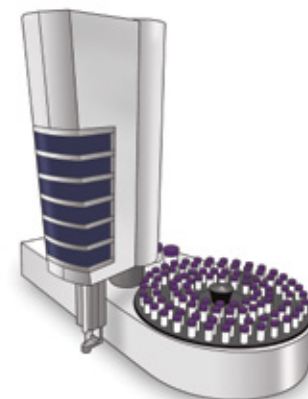

Varian-Chrompack 8400 Autosampler

## Varian-Chrompack Autosampler Syringes

### Varian CP 8400/8410 CP9019/9050 Syringes

Volume	5 µL	10 µL
Model	75	701
Gauge	26s	26s

Cemented Needle (N) Syringe (point style 2)	87900	80300
---	-------	-------

Cemented Needle (N) Syringe 6 pk (point style 2)		80366
--	--	-------

Note: Needles are 51 mm

### Varian 8100/8200 Syringes

Volume	10 µL
Model	701
Gauge	0.48 mm

Removable Needle (RN) Syringe (point style 5)	202880
---	--------

Note: Needles are 50.5 mm





## Shimadzu Autosampler Syringes

### Shimadzu AOC-9 Syringes

Volume	5 $\mu$ L	10 $\mu$ L
Model	75	701
Gauge	26s	26s

Removable Needle (RN) Syringe (point style 2)	87930	80330
--	-------	-------

Note: Needles are 51 mm


### Shimadzu AOC-14/AOC-17/AOC-20 Syringes

Volume	5 $\mu$ L	10 $\mu$ L
Model	75	701
Gauge	22s	22s

Removable Needle (RN) Syringe (point style 2)	202630	202640
--	--------	--------

Cemented Needle (N) Syringe (point style AS)		93898-01
---	--	----------

Note: Needles are 43 mm

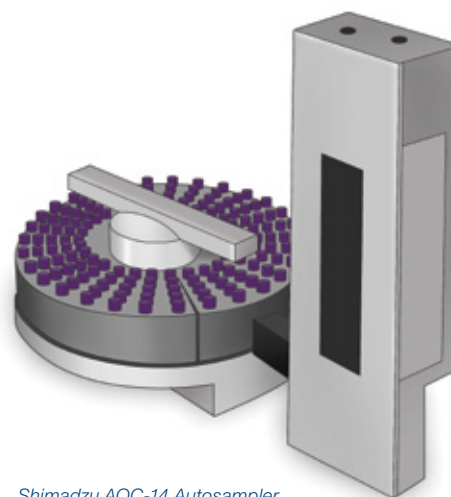

*Shimadzu AOC-14 Autosampler*




# Liquid Chromatography Syringes

Hamilton manufactures a complete range of high performance liquid chromatography (HPLC) syringes and thin layer chromatography (TLC) syringes.

## HPLC Syringes

HPLC syringes are especially designed for manual or automatic injection within high pressure ports. The needle is blunt and electro-polished to slide smoothly into the injection port reducing the wear on critical seals. The selection of the appropriate syringe for HPLC applications takes the following three aspects into account: the type of injection valve, the sample properties and the injection volume. The volume of the syringe should be adjusted to the volume of the injection. In partial loop filling, the injection volume should not be greater than half of the loop capacity. In complete loop filling, where the loop size sets the injection volume, the syringe capacity should be greater than twice the loop volume.

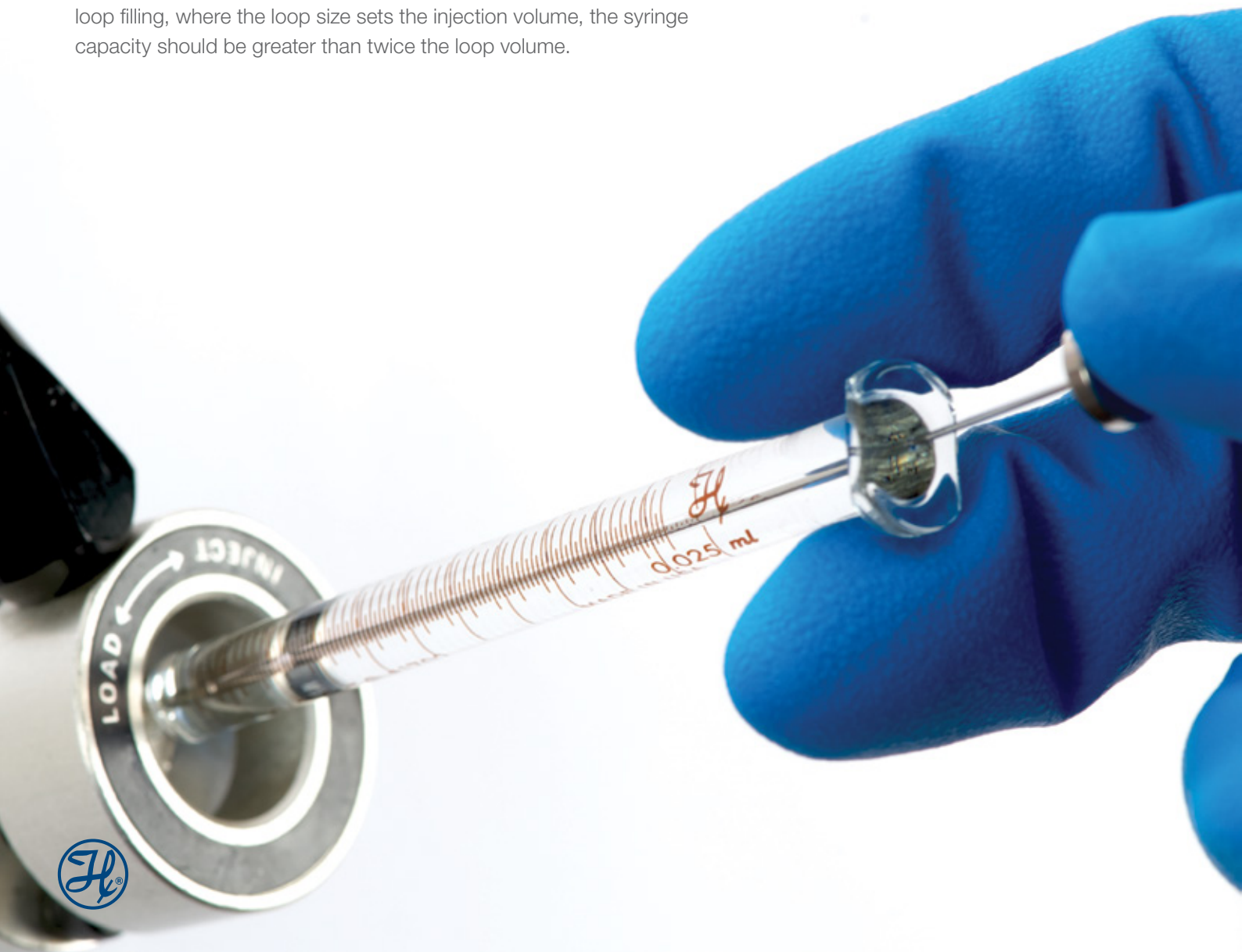

Rheodyne®, Valco®, Beckman Coulter, Altex and SSI Injection Valve Syringes



Gastight Syringes For Aqueous Samples

Volume	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL	1000 µL
Model	1701	1702	1705	1710	1725	1750	1001
Gauge	22s	22s	22s	22s	22s	22	22
Cemented Needle (N) Syringe (point style 3)		80275	80975	81075			
Cemented Needle (LTN) Syringe (point style 3)					81175	81216	81316
Removable Needle (RN) Syringe (point style 3)	80065	80265	80965	81065	81165 <sup>1</sup>	81265	81365
Replacement Needle 6 pk (point style 3)	7770-01	7770-01	7770-01	7770-01	7780-03	7780-04	7780-04

Note: Needles are 51 mm  
 1) 22 gauge needle



Microliter Syringes For Volatile Organic Samples


Volume	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL
Model	701	702	705	710	725	750
Gauge	22s	22s	22s	22s	22	22
Cemented Needle (N) Syringe (point style 3)	80365	80465	80565	80665	80765	80865

Note: Needles are 51 mm



**Tech Tip**

Gastight syringes have a polymer tipped plunger, often PTFE, which essentially wipes the interior of the syringe barrels. This reduces the risk of deposition on the barrel which may cause cross-contamination or plunger seizing.

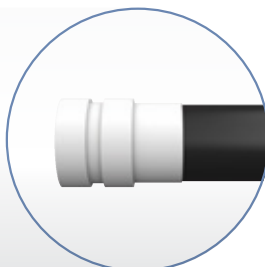

Valco® VISF-1 Valve Syringes



Valco® VISF-1 Injection Valve Syringes

Volume	25 µL	50 µL	100 µL	250 µL	500 µL	1000 µL
Model	1702	1705	1710	1725	1750	1001
Gauge	22s	22s	22s	22s	22	22
Removable Needle (RN) Syringe (point style 3)	80231	80931	81031	81131	81231	81331
Replacement Needle 6 pk (point style 3)	7787-01	7787-01	7787-01	7806-02 <sup>1</sup>	7787-02	7787-02

Note: Needles are 19 mm  
1) Please specify length and point style

Adapter For Repeated Manual Injection

Volume	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL
Model	701/1701	702/1702	705/1705	710/1710	725/1725	750/1750
Reproducibility (Chaney) Adapter	14700	14725	14725	14725	14725	14725



Tech Tip

Improve the reproducibility of multiple manual injections with the Chaney adapter (p/n 14700 and 14725).

# The Power of Polymeric HPLC Columns

Hamilton was the first to offer pressure-stable HPLC columns, and the company now features a wide range of column retention selectivities and performance benefits. Hamilton offers 19 different polymer-based HPLC columns for reversed phase, anion exchange, cation exchange and ion exclusion separations. Hamilton's polymeric HPLC columns combine the inertness and pH stability of polymeric supports with the pressure stability and durability of silica-based columns. Two silica-based C8 and C18 columns are also available for reversed phase separations.



# Thin Layer Chromatography

The point of the needle of a thin layer chromatography (TLC) syringe is covered with a PTFE coating (point style 3T). The coating reduces the surface tension between the needle and the liquid making it ideal for reproducible sample spotting.

## Thin Layer Chromatography Syringes

Volume	10 $\mu$ L	25 $\mu$ L	50 $\mu$ L	100 $\mu$ L
Model	1701	1702	1705	1710
Gauge	26s	22s	22	22
Cemented Needle (N) Syringe (point style 3T)	80050	80250	80950	81050

Note: Needles are 51 mm

For a full list of all syringes available, visit [www.hamiltoncompany.com/syringes](http://www.hamiltoncompany.com/syringes)



## HPLC Autosampler Syringes

HPLC autosamplers enable the automatic introduction of samples into the sample loop as well as a certain amount of sample preparation. Automatic injection has become very common, as it provides improved reproducibility and speed. Hamilton offers a wide range of HPLC autosampler syringes. The features of the autosampler syringes are adapted for an exact fit to a specific autosampler.

### CTC PAL<sup>®</sup> HTC-xt, HTS-xt and HTX-xt Autosampler Syringes

CTC PAL<sup>®</sup> autosamplers are one of the most popular sample processing systems for HPLC analyses. Customers rely on its ability to sample from many different vial and well types and on its ultimate platform flexibility. The syringe is at the heart of every injection that the LC PAL<sup>®</sup> system performs and CTC chose Hamilton to provide this critical component.




## C-Line Syringes

C-Line syringes incorporate several unique design features that ensure superior performance. The innovative direct attachment design of the needle to the barrel eliminates contact between the sample solvent and adhesive reducing carryover and ghost peaks. The inert plunger tip polymer provides enhanced solvent compatibility and longer lifetime.

### C-Line CTC PAL® HTC-xt, HTS-xt and HTX-xt Syringes

Volume	10 µL	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL	5000 µL
Model	1701	1702	1710	1725	1750	1001	1002	1005
Gauge	22s	22s	22s	22	22	22	22	22
Fixed Needle (FN) Syringe (point style 3)	203194 <sup>2</sup>	203075	203077 <sup>2</sup>	203079	203349	203081 <sup>1</sup>	203083 <sup>1</sup>	203085 <sup>1</sup>

Note: Needles are 51 mm  
 1) LTN termination      2) Slimline syringe

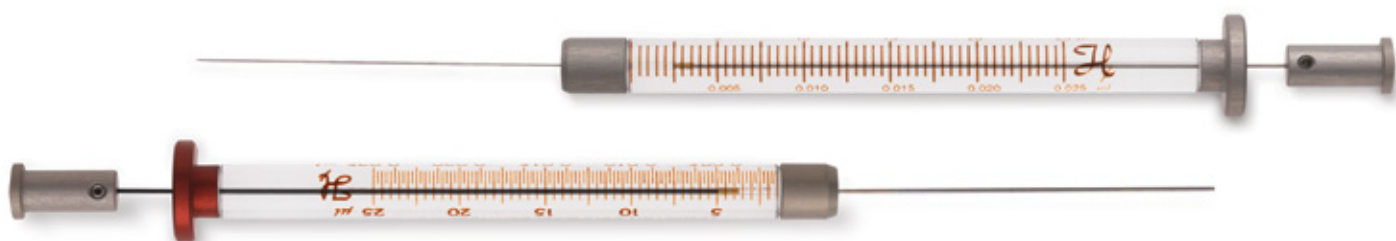

## X-Type Syringes


Hamilton, together with CTC Analytics, has developed X-Type syringes with near zero carryover and a long-lasting plunger tip to meet the requirement of high-throughput life science analytical syringes for sensitive applications. The syringe glass barrel is polished and sealed for inertness and enhanced lifetime. The needle is deactivated to reduce sample adsorption. Tests with critical samples like phospholipids, basic molecules and peptides have shown that carryover and absorption can be reduced depending on the sample. The plunger's life cycle is also improved.

### X-Type Syringes For Sensitive Chromatography Applications

Volume	25 µL	100 µL	100 µL
Model	1702	1710	1710
Gauge	22s	22	22s
Fixed Needle (FN) Syringe (point style 3)	204475	204400	204452

Note: Needles are 51 mm

A guide to syringe terminations and needle points can be found on page 4.

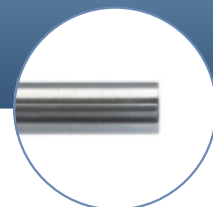

PREP. ROZT. LEU-ENKE  
100 mL A.CN  
100 mL H<sub>2</sub>O  
100 µL HCOOH  
100 µL LEU-ENKE  
(1 mg/mL)





## Agilent/Varian Autosampler Syringes

### Agilent 1290 Infinity LC Injector HTC/HTS Syringes

Volume	10 µL	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL	5000 µL
Model	1701	1702	1710	1725	1750	1001	1002	1005
Gauge	22s	22s	22s	22	22	22	22	22
Fixed Needle (FN) Syringe (point style 3)	203194 <sup>1</sup>	203075	203077 <sup>1</sup>	203079	203349	203081 <sup>2</sup>	203083 <sup>2</sup>	203085 <sup>2</sup>

Note: Needles are 51 mm  
 1) Slimline syringe      2) LTN termination

### Varian Prostar 410/420/430 Syringes

Volume	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL
Model	1702	1710	1725	1750	1001	1002
Gauge	No needle	No needle	No needle	No needle	No needle	No needle
ChemSeal (C) Syringe	62161-01	54658-01	54659-01	54660-01	54661-01	54662-01



Spark Holland Midas Autosampler

## Thermo Scientific® Autosampler Syringes

### AS 100/300 and AS 1000/3000/3500 Syringes

Volume	250 µL	500 µL	1000 µL	2500 µL	5000 µL
Model	1725	1750	1001	1002	1005
Gauge	No needle	No needle	No needle	No needle	No needle
ChemSeal (C) Syringe	202145	202192	81360	81460	81560

### SP 8780/8875/8880 Syringes

Volume	250 µL	2500 µL
Model	1725	1002
Gauge	No needle	No needle
ChemSeal (C) Syringe	202145	81460




Thermo Scientific® AS 3000 Autosampler

### Spark Holland Autosampler Syringes

The Hamilton Spark Holland syringe line works with the following manufacturers' autosamplers:

- ▶ Spark Holland Midas, Triathlon, Endurance and Alias
- ▶ Antec AS 100 and AS 110
- ▶ Dionex ESA Model 542, 540 and 540MT
- ▶ Beckman Coulter 508
- ▶ Grace/Alltech Model 580

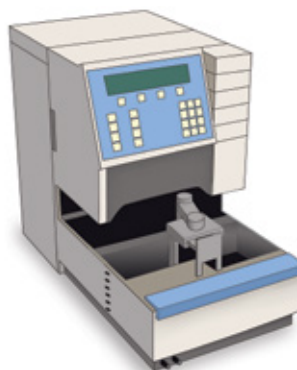


Spark Holland Endurance Autosampler



### Spark Holland Autosampler Syringes

Volume	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL
Model	1702	1710	1725	1750	1001	1002
Gauge	No needle	No needle	No needle	No needle	No needle	No needle
ChemSeal (C) Syringe	62161-01	54658-01	54659-01	54660-01	54661-01	54662-01

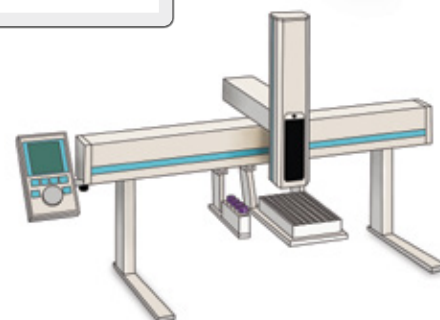

Waters® 2777 Sample Manager

### Waters® Autosampler Syringes

#### Waters® 2777 Sample Manager Syringes

Volume	10 µL	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL	5000 µL
Model	1701	1702	1710	1725	1750	1001	1002	1005
Gauge	22s	22s	22s	22	22	22	22	22
Fixed Needle (FN) Syringe (point style 3)	203194 <sup>1</sup>	203075	203077 <sup>1</sup>	203079	203349	203081 <sup>2</sup>	203083 <sup>2</sup>	203085 <sup>2</sup>

Note: Needles are 51 mm  
 1) Slimline syringes      2) LTN termination

#### Waters® 2707 Autosampler Syringes

Volume	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL
Model	1702	1710	1725	1750	1001	1002
Gauge	No needle	No needle	No needle	No needle	No needle	No needle
ChemSeal (C) Syringe	62161-01	54658-01	54659-01	54660-01	54661-01	54662-01



## PerkinElmer® Autosampler Syringes

## Series 225, Flexar® FX-10/FX-15 Syringes

Volume	25 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL
Model	1702	1710	1725	1750	1001	1002
Gauge	No needle	No needle	No needle	No needle	No needle	No needle
ChemSeal (C) Syringe	62161-01	54658-01	54659-01	54660-01	54661-01	54662-01

## Series 200 Syringes

Volume	50 µL	250 µL	500 µL	1000 µL	2500 µL
Model	1705	1725	1750	1001	1002
Gauge	No needle	No needle	No needle	No needle	No needle
ChemSeal (C) Syringe	80962	81162	81262	81360	81460



PerkinElmer® Series 225 Autosampler



# Sample Preparation and Liquid Handling

Syringe pumps are the most accurate and precise tool for liquid handling. For any syringe pump, the key in achieving the most accurate dispenses is the syringe. Hamilton has a 40-year proven history of designing and manufacturing syringe pumps. With the recent launch of the Microlab® 600, Hamilton's innovation has been brought a step further. The new generation of semi-automated diluters and dispensers are designed for ease of use thanks to the full-color touchscreen interface. The Microlab 600 is used in analytical laboratories performing sample dilution and aliquoting in fields such as pharmaceutical, mining or environmental analysis.



## Syringes For the Hamilton Microlab 600

### Microlab 600 Syringes

Volume	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL
Bubble Free Prime Syringe	59000-05	59000-10	59000-15	59000-20	59000-25	59000-30
Volume	1 mL	2.5 mL	5 mL	10 mL	25 mL	50 mL
Bubble Free Prime Syringe	59000-35	59000-40	59000-45	59000-50	59000-55	59000-60
SaltLine™ Syringe	208335		208336	208337		

### Syringes For Laboratory Syringe Pumps

Laboratory syringe pumps, which may be manufactured by companies such as Harvard Apparatus, KD Scientific, Inc., Kent Scientific Corporation, cetoni and others, are common for microfluidics, drug/nutritional infusions, electrospinning, mass spectrometer calibration, cell feeding and many other applications. There are several infusion pump manufacturers, but one thing is consistent—when a precise flow is critical, all manufacturers turn to Hamilton syringes. Our syringes are available in volumes from 0.5 µL to 100 mL and a variety of terminations to serve almost any application. The PTFE Luer Lock (TLL) termination is extremely common as it provides the most flexible options for connection of tubing.

### Laboratory Syringe Pump Syringes

Volume	50 µL	100 µL	250 µL	500 µL	1 mL	2.5 mL
<b>Model</b>	<b>1705</b>	<b>1710</b>	<b>1725</b>	<b>1750</b>	<b>1001</b>	<b>1002</b>
Luer Lock (TLL) Syringe	80920	81020	81120	81220	81320	81420
SaltLine Luer Lock (TLL) Syringe			203220	203230	203240	203250
Volume	5 mL	10 mL	25 mL	50 mL	100 mL	
<b>Model</b>	<b>1005</b>	<b>1010</b>	<b>1025</b>	<b>1050</b>	<b>1100</b>	
Luer Lock (TLL) Syringe	81520	81620	82520	85020	86020	
SaltLine Luer Lock (TLL) Syringe	203260	203270				

### Tubing Assemblies For Luer Lock and Luer Tip Syringes

Hub Length Material	No hub Custom PTFE	1x Female Luer 300 mm PTFE	1x Female Luer Custom PTFE
Gauge 12	20912	90612	90662
Gauge 18	20918	90618	90668
Gauge 22	20922	90622	90672
Gauge 30	20930	90630	90680



**Tech Tip**

SaltLine

syringes have a highly inert and non-corrosive plunger as well as a particularly resistant ultra-high molecular weight poly-ethylene (UHMWPE) plunger tip, which considerably enhances the syringe lifetime with salt and acid solutions. SaltLine syringes are recommended for life science applications with salty buffers, mining applications with acids or trace metal analysis.

## Automate Your Handheld Pipetting

The Microlab 600 is a highly precise syringe pump that specializes in tasks too small to automate but too large to reliably accomplish by hand, making it invaluable to labs looking to increase efficiency and consistency while reducing cost and wasted buffer. It uses a positive displacement system to provide better than 99% accuracy, independent of a liquid's viscosity, vapor pressure, and temperature—even with harsh chemicals.

With its simple touchscreen interface, built-in wizards and custom method tool, it helps every laboratory work smarter. The free compliance software package on Advanced models features enhanced security options that adhere to CFR 21 Part 11 and FDA GMP/GLP requirements.



For more information on the Microlab 600 benefits your lab, visit [www.hamiltoncompany.com/microlab600](http://www.hamiltoncompany.com/microlab600)



# Life Science Syringes

Hamilton Company offers a wide variety of specialized syringes for life science applications. From custom needles to specialized animal injections and gel loading syringes, Hamilton has a solution for your precision liquid handling life science applications.

## General Animal Injections

Hamilton features a selection of small volume syringes and small gauge needles that can be used for performing injections on lab animals.

### Small Hub Removable Needle Syringes and Custom Needles

Animal injections require a variety of different needle lengths and gauges depending on the target of the injection. For this purpose, Hamilton provides a variety of custom needle options between 27 and 34 gauge. These needles are combined with our small volume Removable Needle syringes to create the perfect application specific syringe. The most common point style for animal injection is point style 4 which comes standard with a sharp 12° beveled tip. If the vessel you are targeting is extremely small, it is possible to select a custom bevel between 12° and 60° to best suit the application. If a blunt point is desired, select point style 3. To minimize needle flexing it is recommended to choose the shortest length possible. The ideal length is typically 13 mm and shorter.

### General Animal Injection Syringes

Volume Model	5 µL	10 µL	25 µL	50 µL	100 µL
Removable Needle (RN) Syringe (no needle included)	7634-01	7653-01	7654-01	7655-01	7656-01

### Custom Stainless Steel Needles

Gauge	27	28	29	30	31	32	33	34
O.D. X I.D. (mm)	0.41 x 0.21	0.36 x 0.18	0.34 x 0.18	0.31 x 0.16	0.26 x 0.13	0.24 x 0.11	0.21 x 0.11	0.16 x 0.05
Small RN Needle	7803-01	7803-02	7803-06	7803-07	7803-03	7803-04	7803-05	207434 <sup>1</sup>


Note: Please specify length and point style when ordering  
<sup>1</sup>) This needle should be used with a Gastight syringe

For a full list of all syringes available, visit [www.hamiltoncompany.com/syringes](http://www.hamiltoncompany.com/syringes)

# Neuroscience Injections

## Introducing Neuros™ Syringes

Hamilton Neuros syringe technology provides unprecedented functionality for controlled animal injections. The Neuros accurately dispenses volumes between 50 nL to 100 µL through an ultrafine needle with a blunt point (point style 3). Developed specifically for neuroscience applications, the Neuros enables the delivery of microvolumes to an exact location while minimizing injection site damage. Neuros syringes come with two types of protective needle sleeves. The sleeve with a blind stop is perfect for cannulated applications and ensures targeted administration with an adjustable penetration depth. The version without a blind stop works best with stereotaxic holders. Both types provide an adjustable needle exposure of 0 to 20 mm.

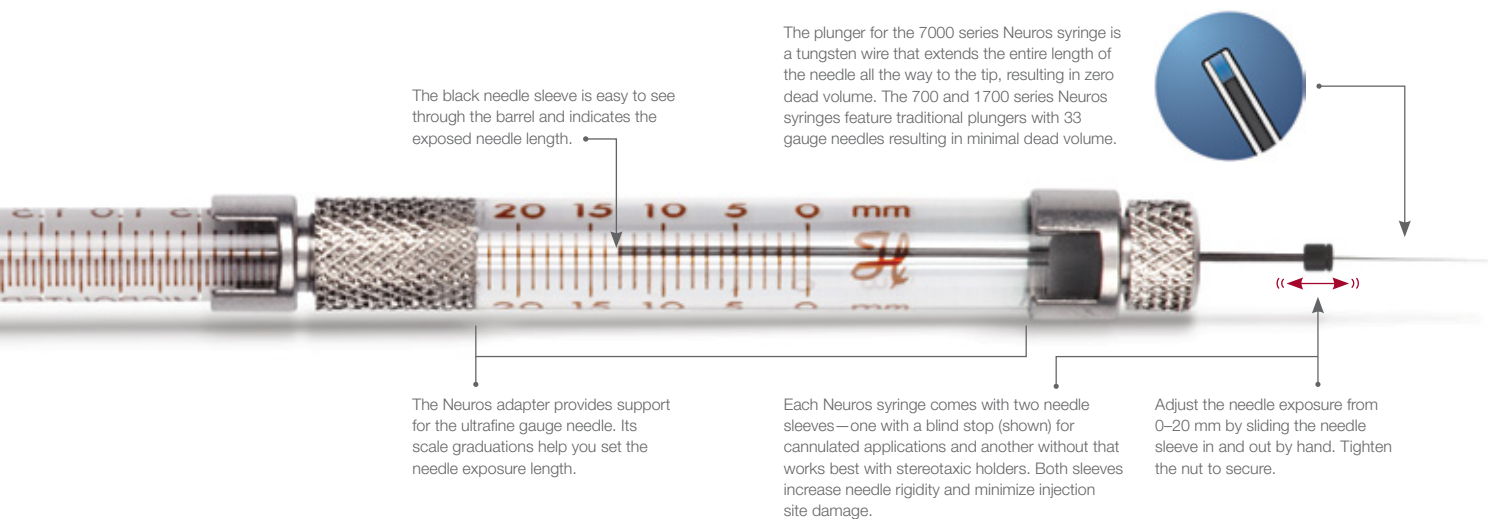

## Neuros Syringes

Volume	0.5 µL	1 µL	2 µL	5 µL	10 µL	25 µL	50 µL	100 µL
Model	7000.5	7001	7002	75	1701	1702	1705	1710
Gauge	32	32	30	33	33	33	33	33
Neuros Syringe (point style 3)	65457-01	65458-01	65459-01	65460-02	65460-05	65460-10	65460-15	65460-20

Note: The exposed needle length is adjustable from 0 to 20 mm

For more information on Neuros syringes, visit [www.hamiltoncompany.com/neuros](http://www.hamiltoncompany.com/neuros)





## Glass Capillary Needles

Glass capillary needles offer several distinct advantages over traditional stainless steel needles. Glass needles can be pulled to form a much finer point, generally smaller than 50  $\mu\text{m}$ . Even with such a small outer diameter, the glass remains rigid and resists dulling. Unlike a beveled needle, the liquid exits out the point of the needle. For several applications these features result in less damage to the animal and an overall better result for the researcher.

To simplify the use of glass capillary needles, Hamilton provides a compression fitting for 1 mm O.D. capillary tubing. When tightened into a small hub Removable Needle syringe, the compression fitting grabs onto the outside of the glass tubing to secure it to the syringe. Next, a priming kit is used to fill the syringe and needle with buffer or mineral oil, and the system is ready for use.

### RN Compression Fitting 1 mm (p/n 55750-01)

Adapts a small hub Removable Needle (RN) termination to a 1 mm pulled glass needle.



## Glass Needle Compression Fittings

Compression Fittings (1 mm Glass Needles)	55750-01
Syringe Priming Kit	PRMKIT

Note: Hamilton does not sell glass needles

### Compatible Syringes For 1 mm Glass Needle Compression Fittings

Volume	5 $\mu\text{L}$	10 $\mu\text{L}$	25 $\mu\text{L}$	50 $\mu\text{L}$	100 $\mu\text{L}$
Model	75	1701	1702	1705	1710
Removable Needle (RN) Syringe (no needle included)	7634-01	7653-01	7654-01	7655-01	7656-01

## Remote Dispensing From a Syringe

There are applications where extending the fluid line between the syringe and needle is required. An example of this is using an infusion pump in combination with a stereotaxic frame. To facilitate this, Hamilton provides compression fittings that are compatible with 1/16 inch PEEK tubing. PEEK tubing is often used for HPLC and was chosen for this purpose because of its high pressure rating. When dispensing extremely small volumes, it is critical to choose tubing that will not flex during the dispense. Flexible tubing has a serious impact on the precision and reproducibility of a dispense.

The compression fitting enables the connection of PEEK tubing to any small hub Removable Needle syringe. Then a small RN coupler connects to the other end of the tubing and to the desired needle. After, the RN coupler can be mounted into the traditional syringe holder for the stereotaxic frame.


### RN Compression Fitting 1/16 inch (p/n 55751-01)

Adapts a small hub Removable Needle (RN) termination to 1/16 inch PEEK tubing.



### Dual Small Hub RN Coupler (p/n 55752-01)

Connects any small hub RN connection to another small hub RN connection.

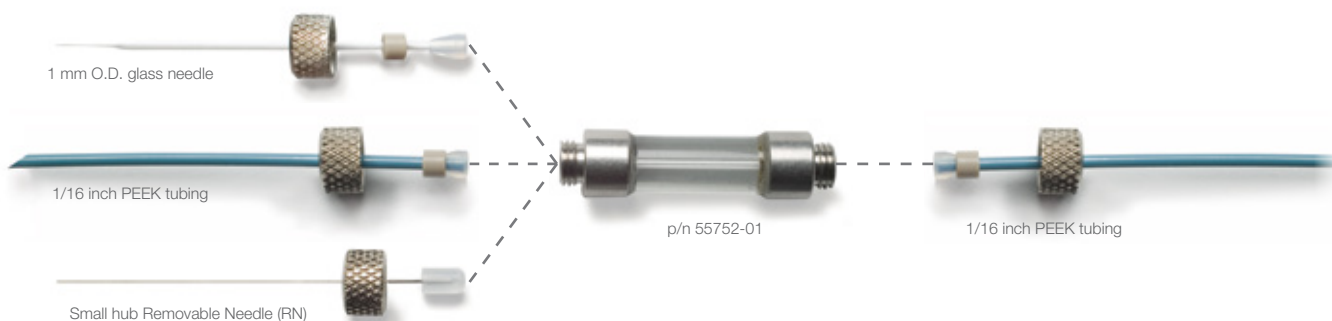

### Compression Fittings For Remote Injections

Compression Fittings (1 mm Glass Needles)	55750-01
Compression Fittings (1/16 inch PEEK Tubing)	55751-01
Dual RN Coupler	55752-01
Syringe Priming Kit	PRMKIT

Note: Hamilton does not sell PEEK tubing or glass needles

### Compatible Syringes For Remote Injection Compression Fittings

Volume	5 $\mu$ L	10 $\mu$ L	25 $\mu$ L	50 $\mu$ L	100 $\mu$ L
Model	75	1701	1702	1705	1710
Removable Needle (RN) Syringe (no needle included)	7634-01	7653-01	7654-01	7655-01	7656-01

### Stereotaxic Infusion Pumps

To improve control over injection rate a variety of stereotaxic mounted infusion pumps are available for Hamilton syringes. These pumps typically fit all standard Hamilton syringes from 0.5 to 500  $\mu$ L and can achieve flow rates as low as 3.66 pL/min. These pumps attach to the flange and barrel allowing for the use of almost any needle length and termination including the new Hamilton Neuros Syringes. Check with the manufacturer of your infusion pump for syringe compatibility.



# Genomics and Biochemistry

## Single and Multi-Channel Gel Loading Syringes

Hamilton Single and Multi-channel Gel Loading Syringes easily and accurately transfer up to eight DNA samples to sequencing gels for electrophoresis, to another microwell plate for duplicate analysis, or to nylon membranes for visualization. A single DNA syringe is also available for individual transfers. Three different needle specifications allow either loading of samples into thin sequencing gel wells that are 0.2, 0.3 or 0.4 mm or larger, or directly onto the top of the sequencing gel. The standard needle spacing of 9 mm allows spotting between the gel plates using most square or sharktooth combs with a microwell format. An adjustable stop collar allows you to pre-set the volume anywhere between 0.2  $\mu\text{L}$  to 10  $\mu\text{L}$ .

### Adjustable Multi-Channel Gel Loading Syringes

Volume	10 $\mu\text{L}$
Gel Thickness	0.2 mm
Well Spacing	9–10.8 mm
.....	
Removable Needle (RN) Eight Channel Syringe (point style 3)	84611
.....	
Replacement Needle 4 pk (point style 3)	78633

### Gel Loading Syringes

Volume	10 $\mu\text{L}$	10 $\mu\text{L}$	10 $\mu\text{L}$
Gel Thickness	0.2 mm	0.3 mm	0.4 mm
Well Spacing	9 mm	9 mm	9 mm
.....			
Cemented Needle (N) Single Channel Syringe (point style 3)	80081	84505	84504
.....			
Cemented Needle (N) Eight Channel Syringe (point style 3)	84511	84503	84502
.....			
Cemented Needle (N) Replacement Syringe	80023	80022	80021



Gel Loading Syringe



## Terasaki Dispenser


The multi-syringe dispenser accurately dispenses six samples simultaneously up to 50 times with the press of a lever. It has been designed for rapid filling of 60 well HLA/Terasaki plates with 6.35 mm well spacing.

### Terasaki Dispenser

Volume	25 µL	50 µL	100 µL	250 µL	500 µL
Model	1702	1705	1710	1725	1750
Volume to Dispense	0.5 µL	1 µL	2 µL	5 µL	10 µL
Removable Needle (RN) Syringe with Dispenser	201900	201910	201920	201930	201940



PB600-1 Repeating Dispenser

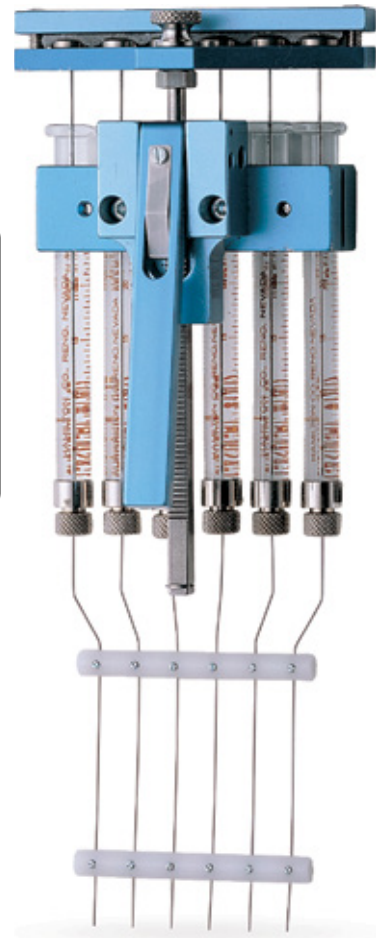

Terasaki Dispenser

## Micro-Pipetting

The PB600-1 Repeating Dispenser consistently dispenses sample between 200 nL–50 µL, up to 50 times with the push of a button. Typical applications include protein crystallization, serology or cell culture.

### Micro-Dispensing

Volume	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL	1000 µL	2500 µL
Model	1701	1702	1705	1710	1725	1750	1001	1002
Gauge	26s	22s	22s	22s	22s	22	22	22
Volume Dispensed Per Click	0.2 µL	0.5 µL	1 µL	2 µL	5 µL	10 µL	20 µL	50 µL
Cemented Needle (N) Syringe (point style 3)	80075	80275	80975	81075	81175 <sup>1</sup>	81216 <sup>1</sup>	81316 <sup>1</sup>	81416 <sup>1</sup>
PB600-1 Repeating Dispenser (no syringe included): 83700 <sup>2</sup>								

Note: Needles are 51 mm

1) LTN termination      2) All syringes listed in this table are compatible with the PB600-1 dispenser

# Product Listing by Features

This section provides feature-based descriptions of Hamilton syringes and needles along with a selection of top items within each category. Product families, series and their components are defined to help you determine the correct item that suits your specific application.


*A guide to syringe terminations and needle points can be found on page 4.*

## Microliter Syringes

Microliter syringes are for use with liquids only. They incorporate a hand-fitted stainless steel plunger with a finely bored syringe barrel. These syringes are ideal for aqueous and organics samples that are not prone to precipitation, crystallization or bonding with glass.

### 700 Series

The stainless steel plunger is manufactured to fit the glass barrel with a tolerance held within millionths of an inch, resulting in unsurpassed syringe life.

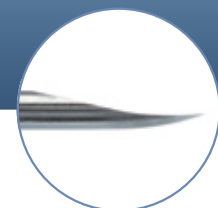

Volume	5 µL	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL
<b>Model</b>	75	701	702	705	710	725	750
<b>Gauge</b>	26s	26s	22s	22s	22s	22s	22
Cemented Needle (N) Syringe (point style 2)	87900	80300 80366 (6 pk)	80400	80500	80600	80700	80800
Cemented Needle (N) Syringe (point style 3)	87919	80383 80365 <sup>2</sup>	80465	80565	80665	80765 <sup>1</sup>	80865
Cemented Needle (N) Syringe (specify needle)	87908	80308	80408	80508	80608	80708	80808
Removable Needle (RN) Syringe (point style 2)	87930	80330 80336 (6 pk)	80430	80530	80630	80730	80830
Removable Needle (RN) Syringe (needle sold separately)	7634-01	7635-01	7636-01	7637-01	7638-01	7639-01	7640-01
Luer Tip (LT) Syringe (needle sold separately)		80301	80401	80501	80601	80701	80801

Note: Needles are 51 mm unless user specified  
 1) 22 gauge                      2) 22s gauge



### 800 Series

The 800 series syringes are a reinforced plunger version of the 700 series. A removable blue syringe holder is attached to the flange of the syringe that prevents heat transfer, dispense inaccuracies, plunger blow outs and plunger bending during injection.




Volume	5 µL	10 µL	25 µL	50 µL	100 µL	250 µL
Model	85	801	802	805	810	825
Gauge	26s	26s	22s	22s	22s	22s
Removable Needle (RN) Syringe (point style 2)	84851	84853	84855	84857	84859	84861
Removable Needle (RN) Syringe (needle sold separately)	7641-01	7642-01	7643-01	7644-01	7645-01	7646-01

Note: Needles are 51 mm

### 7000 Series

The 7000 series syringe is used for dispensing minute liquid volumes and features a Modified Microliter design. The end of the plunger has a tungsten wire extension that fits the entire length of the needle all the way to the tip, resulting in a zero dead volume syringe.

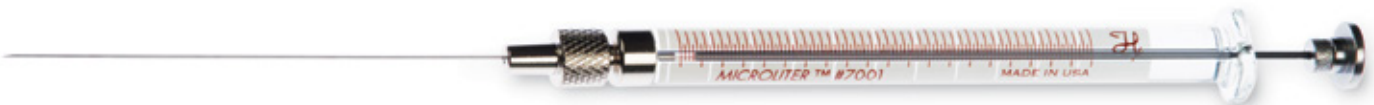

Volume	0.5 µL	0.5 µL	1 µL	1 µL	2 µL	2 µL	5 µL
Model	7000.5	7000.5	7001	7101	7002	7102	7105
Gauge	32	25	25	22	25	23	24
Knurled Hub (KH) Syringe (point style 2)		86259	80135	86211	88411	88511	88011
Knurled Hub (KH) Syringe (point style 3)	86257 <sup>1</sup>	86250	80100	86200	88400	88500	88000


Note: Needles are 70 mm  
1) Needle length is 100 mm

For a full list of all syringes available, visit [www.hamiltoncompany.com/syringes](http://www.hamiltoncompany.com/syringes)

# Gastight Syringes

Gastight syringes are ideal for dispensing both liquids and gases. They have a precision-machined polymer plunger tip, often PTFE, which creates a leak-free seal. With the tight fit, the tip essentially wipes the interior of the syringe barrel free of sample. This feature is useful for samples with dissolved or suspended molecules because it reduces the chance a deposit will occur and result in a damaged plunger.

*A guide to syringe terminations and needle points can be found on page 4.*

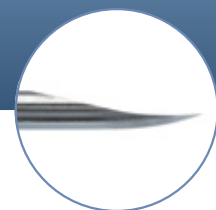

## 1000 Series

The 1000 series syringe is a mid-volume Gastight syringe.



Volume	1 mL	2.5 mL	5 mL	10 mL	25 mL	50 mL	100 mL
<b>Model</b>	<b>1001</b>	<b>1002</b>	<b>1005</b>	<b>1010</b>	<b>1025</b>	<b>1050</b>	<b>1100</b>
<b>Gauge</b>	<b>22</b>	<b>22</b>	<b>22</b>	<b>22</b>	<b>22</b>	<b>22</b>	<b>22</b>
Cemented Needle (LTN) Syringe (point style 2)	81317	81417	81517	81617			
Cemented Needle (LTN) Syringe (point style 3)	81316	81416	81516	81616			
Cemented Needle (LTN) Syringe (specify needle)	81314	81414	81514	81614			
Removable Needle (RN) Syringe (point style 2)	81330 81365 <sup>1</sup>	81430	81530	81630			
Removable Needle (RN) Syringe (needle sold separately)	7649-01	7650-01	7651-01	7652-01			
Luer Tip (LT) Syringe (needle sold separately)	81301	81401	81501	81601			
SampleLock (SL) Syringe (point style 2)	81356	81456	81556	81656	86326	86336	86346
PTFE Luer Lock (TLL) Syringe (needle sold separately)	81320	81420	81520	81620	82520	85020	86020
ChemSeal (C) Syringe (no needle available)	81360	81460	81560	81660			

Note: Needles are 51 mm unless user specified


1) Needle is point style 3





### 1700 Series

The 1700 series syringe is the Gastight version of the 700 series syringes.

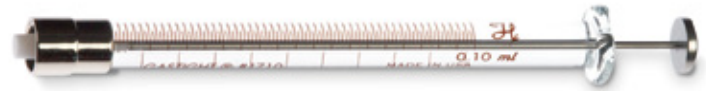

Volume	10 µL	25 µL	50 µL	100 µL	250 µL	500 µL
<b>Model</b>	<b>1701</b>	<b>1702</b>	<b>1705</b>	<b>1710</b>	<b>1725</b>	<b>1750</b>
<b>Gauge</b>	<b>26s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>	<b>22</b>
Cemented Needle (N) Syringe (point style 2)	80000	80200	80900	81000	81100 <sup>1</sup>	81217 <sup>1</sup>
Cemented Needle (N) Syringe (point style 3)	80075	80275	80975	81075	81175 <sup>1</sup>	81216 <sup>1</sup>
Cemented Needle (N) Syringe (specify needle)	80008	80208	80908	81008	81108 <sup>1</sup>	81214 <sup>1</sup>
Removable Needle (RN) Syringe (point style 2)	80030	80230	80930	81030	81130	81230
Removable Needle (RN) Syringe (point style 3)	80065 <sup>2</sup>	80265	80965	81065	81165 <sup>3</sup>	81265
Removable Needle (RN) Syringe (needle sold separately)	7653-01	7654-01	7655-01	7656-01	7657-01	7658-01
Luer Tip (LT) Syringe (needle sold separately)	80001	80201	80901	81001	81101	81201
SampleLock (SL) Syringe (point style 2)			80956	81056	81156	81256
PTFE Luer Lock (TLL) Syringe (needle sold separately)			80920	81020	81120	81220
ChemSeal (C) (no needle available)	80062	80262	80962	81062	81162	81262

Note: Needles are 51 mm unless user specified  
 1) LTN termination    2) 22s gauge    3) 22 gauge

### 1800 Series

The 1800 series syringe is the Gastight version of the 800 series syringes.








Volume	10 µL	25 µL	50 µL	100 µL	250 µL
<b>Model</b>	<b>1801</b>	<b>1802</b>	<b>1805</b>	<b>1810</b>	<b>1825</b>
<b>Gauge</b>	<b>26s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>	<b>22s</b>
Removable Needle (RN) Syringe (point style 2)	84877	84880	84883	84886	84889
Removable Needle (RN) Syringe (needle sold separately)	7659-01	7660-01	7661-01	7662-01	7663-01

Note: Needles are 51 mm unless user specified

# Needles

Hamilton Company offers the widest variety of standard and custom needles in the industry. With needle gauges from 10 to 34 and lengths from 9.53 mm (0.375 inches) to more than 304.79 mm (12 inches), almost any needle is possible. Standard needle length is 51 mm (2 inches). All needles are manufactured from 304 stainless steel.

## Needle Types

Needle Type		Available Gauges	Remarks
Small Hub Removable Needle (RN)		27-34	Small hub removable (RN) needles can be used with removable (RN)-style syringes that are 100 µL and smaller.
		18-26s	
Large Hub Removable Needle (RN)		20-26s	Large hub removable (RN) needles can be used with removable (RN)-style syringes that are 250 µL and larger.
Kel-F Hub Needle (KF)		10-31	Kel-F hub (KF) needles can be used with Luer Tip (LT) and PTFE Luer Lock (TLL) syringes and connectors.
Metal Hub Needle (N)		10-33	Metal hub (N) (Nickel plated brass) needles can be used with TLL syringes and LT or TLL connectors.



Standard Needles


Gauge	O.D. X I.D. (mm)	Point Style	Small RN	Large RN	Metal Hub	Kel-F Hub
33	0.21 x 0.11	3	7762-06 <sup>1</sup>			
32	0.24 x 0.11	3	7762-05			
31	0.26 x 0.13	3	7762-04			
30	0.31 x 0.16	3	7762-03			
28	0.36 x 0.18	3	7762-02			
27	0.41 x 0.21	3	7762-01			
26s	0.47 x 0.13	2	7758-02	7779-02	90039	90139
		3	7768-01	7780-01	91039	90539
		5	7784-07	7784-03	7729-01	7746-12
26	0.46 x 0.26	2	7758-04	7779-04	90026	90126
		3	7768-02	7780-02	91026	90533
		5	7784-08	7784-04	7729-03	7746-10
25	0.52 x 0.26	2			90025	90125
		3			91025	90525
		5			7729-04	7746-09
24	0.57 x 0.31	2			90024	90124
		3			91024	90524
		5			7729-05	7746-08
23	0.64 x 0.34	2			90023	90123
		3			91023	90523
		5			7729-06	7746-07
22s	0.72 x 0.17	2	7758-03	7779-03	90038	90138
		3	7770-01	7780-03	91038	90534
		5	7784-05	7784-01	7729-02	7746-11
22	0.72 x 0.41	2	7758-01	7779-01	90022	90122
		3	7770-02	7780-04	91022	90134
		5	7784-06	7784-02	7729-07	7746-06
21	0.82 x 0.51	2			90021	90121
		3			91021	90521
		5			7729-08	7746-05
20	0.91 x 0.60	2			90020	90120
		3			91020	90520
		5			7729-09	7746-04

Note: All needles come in a six pack and are 51 mm in length  
 1) Length is 38.1 mm


Custom Needles

Gauge	O.D. X I.D. (mm)	Point Style	Small RN	Large RN	Metal Hub	Kel-F Hub
34	0.16 x 0.05	3 or 4	207434			
33	0.21 x 0.11	2, 3 or 4	7803-05		7747-01	
32	0.24 x 0.11	2, 3 or 4	7803-04		7747-02	
31	0.26 x 0.13	2, 3 or 4	7803-03		7748-17	7750-22
30	0.31 x 0.16	2, 3 or 4	7803-07		7748-16	7750-21
29	0.34 x 0.18	2, 3 or 4	7803-06		7748-15	7750-20
28	0.36 x 0.18	2, 3 or 4	7803-02		7748-14	7750-19
27	0.41 x 0.21	2, 3 or 4	7803-01		7748-13	7750-18
26s	0.47 x 0.13	2, 3 or 4	7804-04	7806-04	7748-19	7750-24
		5 or AS	7731-02	7732-04	7751-19	7752-19
26	0.46 x 0.26	2, 3 or 4	7804-03	7806-03	7748-12	7750-17
		5 or AS	7731-01	7732-03	7751-17	7752-17
25s	0.52 x 0.15	2, 3 or 4	7804-10	7806-09	7748-21	7750-26
		5 or AS	7731-06	7732-09	7751-21	7752-21
25	0.52 x 0.26	2, 3 or 4	7804-05	7806-07	7748-11	7750-16
		5 or AS	7731-05	7732-05	7751-16	7752-16
24	0.57 x 0.31	2, 3 or 4	7804-08	7806-06	7748-10	7750-15
		5 or AS	7731-09	7732-07	7751-15	7752-15
23s	0.64 x 0.12	2, 3 or 4	7804-09	7806-08	7748-20	7750-25
		5 or AS	7731-08	7732-08	7751-20	7752-20
23	0.64 x 0.34	2, 3 or 4	7804-07	7806-05	7748-09	7750-14
		5 or AS	7731-07	7732-06	7751-14	7752-14
22s	0.72 x 0.17	2, 3 or 4	7804-02	7806-02	7748-18	7750-23
		5 or AS	7731-04	7732-02	7751-18	7752-18
22	0.72 x 0.41	2, 3 or 4	7804-01	7806-01	7748-08	7750-13
		5 or AS	7731-03	7732-01	7751-13	7752-13
21	0.82 x 0.51	2, 3 or 4	7804-12	7806-11	7748-07	7750-12
		5 or AS	7731-11	7732-11	7751-12	7752-12
20	0.91 x 0.60	2, 3 or 4	7804-11	7806-10	7748-06	7750-11
		5 or AS	7731-10	7732-10	7751-11	7752-11


Note: All needles come in a six pack. Please specify needle length and point style.

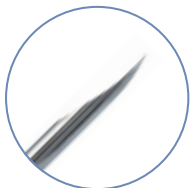

Point style 2

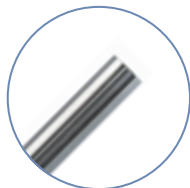


Point style 3

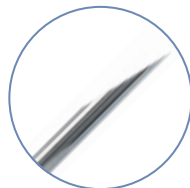

Point style 4



Point style 5

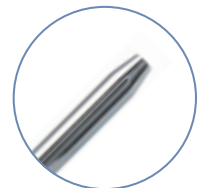

Point style AS



Custom Needles continued

Gauge	O.D. X I.D. (mm)	Point Style	Small RN	Metal Hub	Kel-F Hub
19	1.07 x 0.69	2, 3 or 4	207419	7748-05	7750-10
		5 or AS		7751-10	7752-10
18	1.27 x 0.84	2, 3 or 4	7804-06	7748-04	7750-09
		5 or AS		7751-09	7752-09
17	1.47 x 1.07	2, 3 or 4		7748-03	7750-08
		5 or AS		7751-08	7752-08
16	1.65 x 1.19	2, 3 or 4		7748-02	7750-07
		5 or AS		7751-07	7752-07
15	1.83 x 1.37	2, 3 or 4		7748-01	7750-06
		5 or AS		7751-06	7752-06
14	2.11 x 1.60	2, 3 or 4		7749-05	7750-05
		5 or AS		7751-05	7752-05
13	2.41 x 1.80	2, 3 or 4		7749-04	7750-04
		5 or AS		7751-04	7752-04
12	2.77 x 2.16	2, 3 or 4		7749-03	7750-03
		5 or AS		7751-03	7752-03
11	3.05 x 2.39	2, 3 or 4		7749-02	7750-02
		5 or AS		7751-02	7752-02
10	3.40 x 2.70	2, 3 or 4		7749-01	7750-01
		5 or AS		7751-01	7752-01

Note: All needles come in a six pack. Please specify needle length and point style.

Custom needles usually ship within 5-7 business days.



# Syringe Care and Use

The life of your Hamilton syringe is directly related to its cleanliness and proper care. In general, solvents suitable for routine cleaning include methanol, acetonitrile and acetone. Use solvents of high purity grade. Halogenated hydrocarbons should not be used because they may damage some glue joints.

## Chemical Compatibility

The adhesive used to affix needles and hubs to Hamilton Gastight and Microliter syringes is the most chemically resistant available. However, with prolonged exposure, some solvents may attack and deteriorate this highly resistant adhesive. In particular, caution should be exercised with solvents containing halogenated hydrocarbons such as dichloromethane (methylene chloride). For applications using these solvents, Removable Needle (RN) syringes are recommended because no adhesive is present in the fluid path. Be sure to rinse the syringe thoroughly after each use with a solvent that is known to solubilize your sample followed by a solvent such as acetone to ensure that the glue does not remain in contact with a potentially harmful solvent.

## Cleaning Syringes

Rinse the syringe after use with an appropriate solvent (i.e., a solvent known to be effective in solvating the sample). Preferred cleaning agents are non-alkaline, non-phosphate and non-detergent based. A biodegradable, non-phosphate, organic Cleaning Solution Concentrate is available from Hamilton (p/n 18311). Following the use of a cleaning agent, rinse the syringe with deionized water and finally acetone. Wipe the exterior surfaces of the syringe barrel and needle dry with a lint-free tissue. Make sure that there is no residual cleaning agent in the syringe before using or storing the syringe. Do not soak or submerge the entire syringe in any cleaning agent. Prolonged contact with cleaning agents may damage bonded parts. It is recommended that the syringe be placed into the original packaging after it is properly cleaned.



*Needle Cleaning Kit*



## Clogged Needles

Do not attempt to clean clogged needles by forcing liquid or compressed air through the syringe. Excessive pressure will split the glass barrel. Instead, use the Hamilton Needle Cleaning Kit (p/n 76620). Start by using the cleaning wires to dislodge any foreign material. Then flush with the Cleaning Solution Concentrate to further dissolve the clog. Once the clog is removed, rinse the syringe and needle thoroughly with deionized water. Wipe the exterior surfaces of the syringe barrel and needle dry with a lint-free tissue. Make sure that there is no residual cleaning agent in the syringe before using or storing the syringe.

## Storage

After Hamilton syringes are cleaned and dried store them in the original packaging or on the syringe rack (p/n 204880) to prevent breakage or damage.

## Sterilizing Syringes

Chemical sterilizers are the best means for sterilizing syringes. Hamilton syringes may be sterilized with ethylene oxide. LT, TLL and disassembled RN syringes may be autoclaved at a temperature of up to 115 °C. However, repeated autoclaving will shorten syringe life. When autoclaving is required, remove the plunger from the syringe. Autoclaving cemented needle syringes is not recommended as glass and metal expands at different rates which would compromise the glued connection. If your application only requires disinfection, then Hamilton recommends the use of Microside SQ®.

## Lubricating Syringes

A clean syringe does not require any lubricating grease. Grease may lead to a variety of problems including sample contamination, seizing of the plunger and barrel damage.



Syringe Rack



**Syringe Troubleshooting**

Problem	Possible Cause	Suggested Solution
Clogged needle	Worn GC septa	Replace GC septa
	Poor syringe cleaning	Ensure proper needle cleaning with Hamilton cleaning wires and cleaning concentrate
	Viscous sample	Select a larger needle gauge for viscous sample
Bent plunger	Inappropriate handling technique	Use syringe with reinforced plunger
	Poor storage	Store syringes on the syringe rack or in their original packaging
Seized Microliter plunger	Deposits between the glass barrel and the plunger wire	Avoid operating the syringe dry
		Ensure proper syringe cleaning, it may be possible to free the plunger using a solvent known to solubilize your sample
Low reproducibility	Inappropriate handling technique	Use Reproducibility (Chaney) Adapter
Air bubbles	Insufficient priming	Fully prime the syringe by inserting the needle into the sample and slowly aspirating and dispensing
Needle leakage	Inappropriate sealing of the needle	Remount needle checking that the cleaning wire has been removed and there is no damage to the plastic seal
Cemented needle is loose	Solvents have attacked the glue holding the needle	Use RN syringes that have no glue in the fluid path, especially with halogenated solvents
Gastight plunger leakage	Plunger tip damaged due to age	Replace plunger
	Plunger tip damage due to deposits on the glass	Improve cleaning to eliminate deposits. For some applications SaltLine syringes may provide longer life.
Carryover	Inappropriate syringe cleaning	Ensure proper syringe and needle cleaning
	Inappropriate syringe selection	Use a Gastight syringe with aqueous samples





# Part Number Index

PRMKIT.....	31, 33	7746-11.....	41	7751-12.....	42	7803-06.....	29, 42
7634-01.....	29, 31, 33, 36	7746-12.....	41	7751-13.....	42	7803-07.....	29, 42
7635-01.....	36	7747-01.....	42	7751-14.....	42	7804-01.....	42
7636-01.....	36	7747-02.....	42	7751-15.....	42	7804-02.....	42
7637-01.....	36	7748-01.....	43	7751-16.....	42	7804-03.....	42
7638-01.....	36	7748-02.....	43	7751-17.....	42	7804-04.....	42
7639-01.....	36	7748-03.....	43	7751-18.....	42	7804-05.....	42
7640-01.....	36	7748-04.....	43	7751-19.....	42	7804-06.....	43
7641-01.....	37	7748-05.....	41, 43	7751-20.....	42	7804-07.....	42
7642-01.....	37	7748-06.....	42	7751-21.....	42	7804-08.....	42
7643-01.....	37	7748-07.....	42	7752-01.....	43	7804-09.....	42
7644-01.....	37	7748-08.....	42	7752-02.....	43	7804-10.....	42
7645-01.....	37	7748-09.....	42	7752-03.....	43	7804-11.....	42
7646-01.....	37	7748-10.....	42	7752-04.....	43	7804-12.....	42
7649-01.....	38	7748-11.....	42	7752-05.....	43	7806-01.....	42
7650-01.....	38	7748-12.....	42	7752-06.....	43	7806-02.....	18, 42
7651-01.....	38	7748-13.....	42	7752-07.....	43	7806-03.....	42
7652-01.....	38	7748-14.....	42	7752-08.....	43	7806-04.....	42
7653-01.....	29, 31, 33, 39	7748-15.....	42	7752-09.....	43	7806-05.....	42
7654-01.....	29, 31, 33, 39	7748-16.....	42	7752-10.....	43	7806-06.....	42
7655-01.....	29, 31, 33, 39	7748-17.....	42	7752-11.....	42	7806-07.....	42
7656-01.....	29, 31, 33, 39	7748-18.....	42	7752-12.....	42	7806-08.....	42
7657-01.....	39	7748-19.....	42	7752-13.....	42	7806-09.....	42
7658-01.....	39	7748-20.....	42	7752-14.....	42	7806-10.....	42
7659-01.....	39	7748-21.....	42	7752-15.....	42	7806-11.....	42
7660-01.....	39	7749-01.....	43	7752-16.....	42	14700.....	9, 18
7661-01.....	39	7749-02.....	43	7752-17.....	42	14725.....	9, 18
7662-01.....	39	7749-03.....	43	7752-18.....	42	18311.....	44
7663-01.....	39	7749-04.....	43	7752-19.....	42	20912.....	27
7729-01.....	41	7749-05.....	43	7752-20.....	42	20918.....	27
7729-02.....	41	7750-01.....	43	7752-21.....	42	20922.....	27
7729-03.....	41	7750-02.....	43	7758-01.....	41	20930.....	27
7729-04.....	41	7750-03.....	43	7758-02.....	7, 8, 41	54658-01.....	23, 24, 25
7729-05.....	41	7750-04.....	43	7758-03.....	7, 8, 41	54659-01.....	23, 24, 25
7729-06.....	41	7750-05.....	43	7758-04.....	41	54660-01.....	23, 24, 25
7729-07.....	41	7750-06.....	43	7762-01.....	41	54661-01.....	23, 24, 25
7729-08.....	41	7750-07.....	43	7762-02.....	41	54662-01.....	23, 24, 25
7729-09.....	41	7750-08.....	43	7762-03.....	41	55750-01.....	31, 33
7731-01.....	42	7750-09.....	43	7762-04.....	41	55751-01.....	33
7731-02.....	42	7750-10.....	43	7762-05.....	41	55752-01.....	33
7731-03.....	42	7750-11.....	42	7762-06.....	41	59000-05.....	27
7731-04.....	42	7750-12.....	42	7768-01.....	41	59000-10.....	27
7731-05.....	42	7750-13.....	42	7768-02.....	41	59000-15.....	27
7731-06.....	42	7750-14.....	42	7770-01.....	17, 41	59000-20.....	27
7731-07.....	42	7750-15.....	42	7770-02.....	41	59000-25.....	27
7731-08.....	42	7750-16.....	42	7779-01.....	7, 8, 41	59000-30.....	27
7731-09.....	42	7750-17.....	42	7779-02.....	41	59000-35.....	27
7731-10.....	42	7750-18.....	42	7779-03.....	7, 8, 41	59000-40.....	27
7731-11.....	42	7750-19.....	42	7779-04.....	41	59000-45.....	27
7732-01.....	42	7750-20.....	42	7780-01.....	41	59000-50.....	27
7732-02.....	42	7750-21.....	42	7780-02.....	41	59000-55.....	27
7732-03.....	42	7750-22.....	42	7780-03.....	17, 41	59000-60.....	27
7732-04.....	42	7750-23.....	42	7780-04.....	17, 41	62161-01.....	23, 24, 25
7732-05.....	42	7750-24.....	42	7784-01.....	41	65457-01.....	30
7732-06.....	42	7750-25.....	42	7784-02.....	41	65458-01.....	30
7732-07.....	42	7750-26.....	42	7784-03.....	41	65459-01.....	30
7732-08.....	42	7751-01.....	43	7784-04.....	41	65460-02.....	30
7732-09.....	42	7751-02.....	43	7784-06.....	41	65460-05.....	30
7732-10.....	42	7751-03.....	43	7784-07.....	41	65460-10.....	30
7732-11.....	42	7751-04.....	43	7784-08.....	41	65460-15.....	30
7746-04.....	41	7751-05.....	43	7787-01.....	18	65460-20.....	30
7746-05.....	41	7751-06.....	43	7787-02.....	18	76620.....	45
7746-06.....	41	7751-07.....	43	7803-01.....	29, 42	78633.....	34
7746-07.....	41	7751-08.....	43	7803-02.....	29, 42	80000.....	8, 39
7746-08.....	41	7751-09.....	43	7803-03.....	29, 42	80001.....	39
7746-09.....	41	7751-10.....	43	7803-04.....	29, 42	80008.....	39
7746-10.....	41	7751-11.....	42	7803-05.....	29, 42	80021.....	34

80022	34	80956	39	84503	34	90534	41
80023	34	80962	25, 39	84504	34	90539	41
80030	8, 39	80965	17, 39	84505	34	90612	27
80050	19	80975	17, 35, 39	84511	34	90618	27
80062	39	81000	8, 39	84611	34	90622	27
80065	17, 39	81001	39	84851	37	90630	27
80075	35, 39	81008	39	84853	37	90662	27
80079	13	81020	27, 39	84855	37	90668	27
80080	13	81030	8, 39	84857	37	90672	27
80081	34	81031	18	84859	37	90680	27
80094	13	81050	19	84861	37	91020	41
80096	13	81056	39	84877	8, 39	91021	41
80100	37	81062	39	84880	8, 39	91022	41
80135	7, 37	81065	17, 39	84883	8, 39	91023	41
80176	13	81075	17, 35, 39	84886	8, 39	91024	41
80200	8, 39	81100	8, 39	84889	8, 39	91025	41
80201	39	81101	39	85020	27, 38	91026	41
80208	39	81108	39	86020	27, 38	91038	41
80230	8, 39	81120	27, 39	86200	37	91039	41
80231	18	81130	8, 39	86201	7	201900	35
80250	19	81131	18	86211	37	201910	35
80262	39	81156	39	86250	37	201920	35
80265	17, 39	81162	25, 39	86257	37	201930	35
80275	17, 35, 39	81165	17, 39	86259	7, 37	201940	35
80300	7, 14, 36	81175	17, 35, 39	86274	13	202145	23
80301	36	81201	39	86276	13	202192	23
80308	10, 36	81214	39	86326	38	202630	15
80318	14	81216	17, 35, 39	86336	38	202640	15
80330	7, 15, 36	81217	8, 39	86346	38	202660	14
80336	7, 36	81220	27, 39	87402	10	202880	14
80365	17, 36	81230	8, 39	87900	7, 14, 36	203043	12
80366	7, 14, 36	81231	18	87908	36	203075	21, 23, 24
80383	36	81256	39	87919	36	203076	12
80387	13	81262	25, 39	87930	7, 15, 36	203077	21, 23, 24
80388	13	81265	17, 39	87987	13	203078	12
80389	13	81301	38	87988	13	203079	21, 23, 24
80390	13	81314	38	87989	13	203080	12
80391	13	81316	17, 35, 38	87990	13	203081	21, 23, 24
80393	13	81317	8, 38	87993	13	203082	12
80400	7, 36	81320	27, 38	87994	13	203083	21, 23, 24
80401	36	81330	8, 38	88000	37	203084	12
80408	36	81331	18	88011	7, 37	203085	21, 23, 24
80430	7, 36	81356	38	88400	37	203086	12
80465	17, 36	81360	23, 25, 38	88411	7, 37	203141	12
80500	7, 36	81365	17, 38	88500	37	203181	12
80501	36	81401	38	88511	37	203182	12
80508	36	81414	38	90020	41	203185	12
80530	7, 36	81416	35, 38	90021	41	203189	12
80565	17, 36	81417	8, 38	90022	41	203194	21, 23, 24
80600	7, 36	81420	27, 38	90023	41	203205	12
80601	36	81430	8, 38	90024	41	203220	27
80608	36	81456	38	90025	41	203230	27
80630	7, 36	81460	23, 25, 38	90026	41	203240	27
80665	17, 36	81501	38	90038	41	203250	27
80700	7, 36	81514	38	90039	41	203260	27
80701	36	81516	38	90120	41	203270	27
80708	36	81517	8, 38	90121	41	203349	21, 23, 24
80730	7, 36	81520	27, 38	90122	41	203361	12
80765	17, 36	81530	8, 38	90123	41	204000	14
80800	7, 36	81556	38	90124	41	204001	14
80801	36	81560	23, 38	90125	41	204051	14
80808	36	81601	38	90126	41	204052	14
80830	7, 36	81614	38	90134	41	204400	21
80865	17, 36	81616	38	90138	41	204452	21
80900	8, 39	81617	8, 38	90139	41	204475	21
80901	39	81620	27, 38	90520	41	204880	45
80908	39	81630	8, 38	90521	41	207419	43
80920	27, 39	81656	38	90523	41	207434	29, 42
80930	8, 39	81660	38	90524	41	208335	27
80931	18	82520	27, 38	90525	41	208336	27
80950	19	84502	34	90533	41	208337	27



# About Hamilton Company


Hamilton Company is a global enterprise with affiliates in Reno, Nevada; Franklin, Massachusetts; and Bonaduz, Switzerland and sales offices throughout the world.

We are an industry leader in the design and manufacture of liquid handling, process analytics, robotics and automated storage solutions. For more than 60 years, Hamilton has been satisfying customer needs by combining quality materials with skilled workmanship to ensure the highest level of performance. Hamilton's lifelong commitment to precision and quality has earned us global ISO 9001 Certification.



Founded on the technology of analytical Microliter™ and Gastight® syringes, Hamilton Company has a broad offering of laboratory products including manual and semi-automated precision fluid measuring instruments, chromatography products, process sensors, laboratory electrodes, pipettes and more. Top innovations from these lines include Arc® pH, DO and Conductivity Intelligent Sensors, Microlab® 600 Diluters/Dispensers and the Microlab® 300 Guided Pipetting System.

A pioneer in liquid handling equipment and laboratory automation technology, Hamilton Robotics is known for advancing life science and biotechnology industries through reliability, performance and flexibility. Hamilton Robotics is the industry leader in design and manufacturing with patented technologies such as Compression-induced O-Ring Expansion (CO-RE®), Total Aspiration and Dispensing Monitoring (TADM®) and Anti-Droplet Control (ADC™). Hamilton Robotics' platforms include Hamilton VANTAGE™ Liquid Handling System, its newest vertically-integrated liquid handler, Microlab STAR™, Hamilton Robotics' highest selling automated pipetting platform, and Microlab NIMBUS®, the first in its class of compact, high-speed, personalized pipetting workstations.

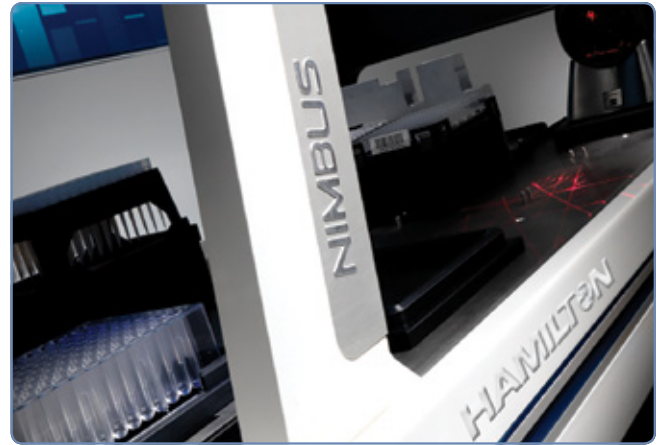

Hamilton Storage Technologies offers comprehensive ultra-low temperature automated sample management systems for microtube and microplate storage. Hamilton Storage Technologies' line of biobanking and compound storage solutions, as well as consumables, are designed for a broad array of life science processes. Products include Hamilton BiOS®, SAM™ and ASM™, designed for sample integrity, flexibility and reliability.

Hamilton Company is focused on blending invention and accuracy to deliver customers unparalleled products.



Your Hamilton Representative

 **ENRICO BRUNO s.r.l.**  
APPARECCHI SCIENTIFICI  
SANITARI – INDUSTRIALI

© 2015 Hamilton Company. All rights reserved.  
Thermo Scientific is a registered trademark of Thermo Fisher Scientific Inc.  
Valco® is a registered trademark of Valco Instruments Co. Inc.  
Rheodyne is a registered trademark of Units of IDEX® Corporation.  
Waters is a registered trademark of Waters.  
PerkinElmer and Flexar are registered trademarks of PerkinElmer, Inc.  
Merlin MicroSeal is a trademark of Merlin Instrument Company.  
Microside SQ is a registered trademark of Global Biotechnologies, Inc. in the U.S. and/or other countries.  
PAL is a registered trademark of CTC Analytics AG, Switzerland.

All other trademarks are owned and/or registered by Hamilton Company in the U.S. and/or other countries.  
Lit. No. L20078 Rev. B – 1/2015