



The price leader in analytical balances with internal or external calibration – now as version with [Max] 220 g!

Features

- New model with incredibly high resolution. Ideal if extremely precise weighing is required
- KERN ADJ: Automatic internal adjustment in the case of a change in temperature $\geq 2\text{ }^{\circ}\text{C}$ or timecontrolled every 3 h, guarantees high degree of accuracy and makes the balance independent of its location of use
- KERN ADB: Adjusting program CAL for quick setting of the balance accuracy using an external *test weight*
- **1** ADB 600-C3: Compact, space-saving carat balance with a readout of 0.001 ct and a weighing range of 600 ct. The high level of accuracy saves hard cash wherever you are weighing valuable precious stones
- Level indicator and levelling feet for precise levelling of the scale, fitted as standard, to give the most accurate weighing result
- Large glass draught shield with 3 sliding doors for easy access to the items being weighed
- Compact size, practical for small spaces
- Simple and convenient 6-key operation

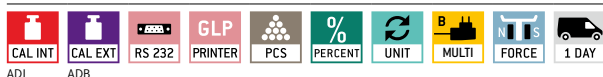
Technical data

- Large backlit LCD display, digit height 16 mm
- Dimensions weighing surface, stainless steel, $\varnothing 90\text{ mm}$
- Overall dimensions (incl. draught shield) W×D×H
KERN ADB/ADJ: 230×310×330 mm
KERN ADB-C/ADJ-C: 230×310×210 mm
- Weighing space W×D×H
KERN ADB/ADJ: 170×160×205 mm
KERN ADB-C/ADJ-C: 170×160×110 mm
- Permissible ambient temperature 10 °C/30 °C

Accessories

- **2** Set for density determination of liquids and solids with density $\leq/\geq 1$, the density is indicated directly on the display, KERN YDB-03
- **3** Ioniser to neutralise electrostatic charge, KERN YBI-01A
- **4** Gemstones plate, aluminium with practical spout, W×D×H 83×66×23 mm, KERN AEJ-A05
- **5** Weighing table to absorb vibrations and oscillations, which would otherwise distort the weighing result, KERN YPS-03
- Minimum weight of sample, smallest weight to be weighed, depending on the required process accuracy, only in combination with a DAkkS calibration certificate, KERN 969-103
- Equipment qualification: compliant qualification concept which includes the following validation services: Installation Qualification (IQ), Operating Qualification (OQ), Further details see 208
- Further details, plenty of further accessories and suitable printers see *Accessories*

STANDARD



OPTION
























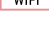


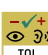




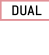




Model	Weighing capacity [Max] g	Readability [d] mg	Reproducibility mg	Linearity mg	Net weight approx. kg	Option	
						DAkkS Calibr. Certificate	
KERN						DAkkS KERN	
ADB 100-4	120	0,1	0,2	$\pm 0,4$	4,4	963-101	
ADB 200-4	220	0,1	0,2	$\pm 0,4$	4,4	963-101	
ADB 600-C3	600 ct	0,001 ct	0,002 ct	$\pm 0,004\text{ ct}$	3,8	963-101	
ADJ							
ADJ 100-4	120	0,1	0,2	$\pm 0,4$	5	963-101	
ADJ 200-4	210	0,1	0,2	$\pm 0,4$	5	963-101	
ADJ 600-C3	600 ct	0,001 ct	0,002 ct	$\pm 0,004\text{ ct}$	4,5	963-101	

New model

= Carat balances

Pictograms

 Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)	 KERN Communication Protocol (KCP): It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems	 Suspended weighing: Load support with hook on the underside of the balance
 Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required		 Battery operation: Ready for battery operation. The battery type is specified for each device
 Easy Touch: Suitable for the connection, data transmission and control through PC or tablet.	 GLP/ISO log: The balance displays serial number, user ID, weight, date and time, regardless of a printer connection	 Rechargeable battery pack: Rechargeable set
 Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	 GLP/ISO log: With weight, date and time. Only with KERN printers.	 Universal plug-in power supply: with universal input and optional input socket adapters for A) EU, CH, GB; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS
 Alibi memory: Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.	 Piece counting: Reference quantities selectable. Display can be switched from piece to weight	 Plug-in power supply: 230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available
 Data interface RS-232: To connect the balance to a printer, PC or network	 Recipe level A: The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out	 Integrated power supply unit: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
 RS-485 data interface: To connect the balance to a printer, PC or other peripherals. Suitable for datatransfer over large distances. Network in bus topology is possible	 Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display	 Weighing principle: Strain gauges: Electrical resistor on an elastic deforming body
 USB data interface: To connect the balance to a printer, PC or other peripherals	 Totalising level A: The weights of similar items can be added together and the total can be printed out	 Weighing principle: Tuning fork: A resonating body is electromagnetically excited, causing it to oscillate
 Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals	 Percentage determination: Determining the deviation in % from the target value (100 %)	 Weighing principle: Electromagnetic force compensation: Coil inside a permanent magnet. For the most accurate weighings
 WiFi data interface: To transfer data from the balance to a printer, PC or other peripherals	 Weighing units: Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details	 Weighing principle: Single cell technology: Advanced version of the force compensation principle with the highest level of precision
 Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.	 Weighing with tolerance range: (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model	 Verification possible: The time required for verification is specified in the pictogram
 Analogue interface: to connect a suitable peripheral device for analogue processing of the measurements	 Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value	 DAKkS calibration possible (DKD): The time required for DAKkS calibration is shown in days in the pictogram
 Interface for second balance: For direct connection of a second balance	 Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.	 Factory calibration (ISO): The time required for Factory calibration is shown in days in the pictogram
 Network interface: For connecting the scale to an Ethernet network		 Package shipment: The time required for internal shipping preparations is shown in days in the pictogram
		 Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg – 2500 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg – 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages DE, EN, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

Your KERN specialist dealer:

 **ENRICO BRUNO s.r.l.**
APPARECCHI SCIENTIFICI
SANITARI – INDUSTRIALI

www.enrico-bruno.it