BRAND Inserts

Genuine innovation for cell culture



High flexibility – to be used as a hanging or standing insert

cellGrade™ plus treated membranes for improved cell attachment and growth

www.enrico-bruno.it

Easy to use and efficient in your cell culture lab



BRAND Inserts - genuine innovation

■ Flexible

Keep using the multiwell plates or dishes you already use in your cell culture lab – no change of cell culture surface when switching from hanging to standing inserts!


When used as a hanging insert, the BRAND Inserts fit nearly all multiwell plates. For bigger wells and culture dishes, simply use the same insert as a standing insert.

■ Functional

The cellGrade™ plus surface of the membranes ensures excellent cell attachment and growth. The highly transparent membranes give you good visibilty of your cell cultures, and the high pore density contributes to good permeability.


■ Easy to use and efficient

Innovative packaging provides for efficient unpacking and reduces waste.

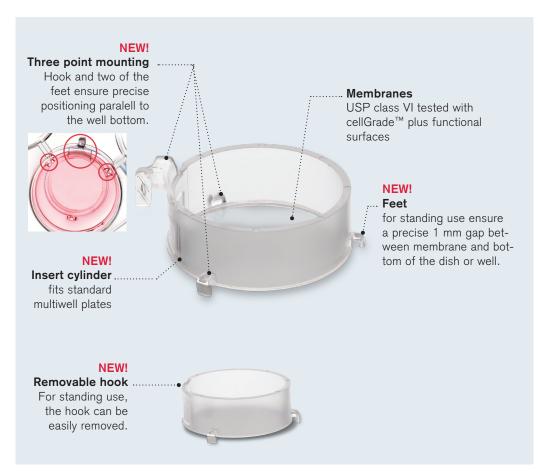

Easy-to-use packaging:

(1) Simply open the entire blister by pulling the sealing paper in one diagonal move (lilac arrow).

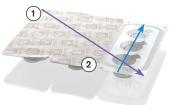

2 A perforation in the sealing paper also allows to open only one line of inserts and leave the others securely sealed (blue arrow).

BRAND Cell culture inserts are available for 6-, 12- and 24-well sizes and in two pore sizes (0.4 μm and 8.0 $\mu m)$ with PC-membrane or PET-membrane.

Individually packed in single blisters or in a plate-style blister.





Distribuito da:



